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Aurum ProTM

Familiar Sequence, Different Flexibility

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UNMATCHED FLEXIBILITY
Heat-treated NiTi gives amazing flexibility, inhibiting transportation of root canal


THE SYSTEM YOU TRUST
The most commonly used rotary file system

INCREDIBLE STRENGTH
20% Resistant to Cyclic Fatigue than Protaper Universal

Amazing Strength of Aurum Pro

Number of cycle for fracture, Internal testing data

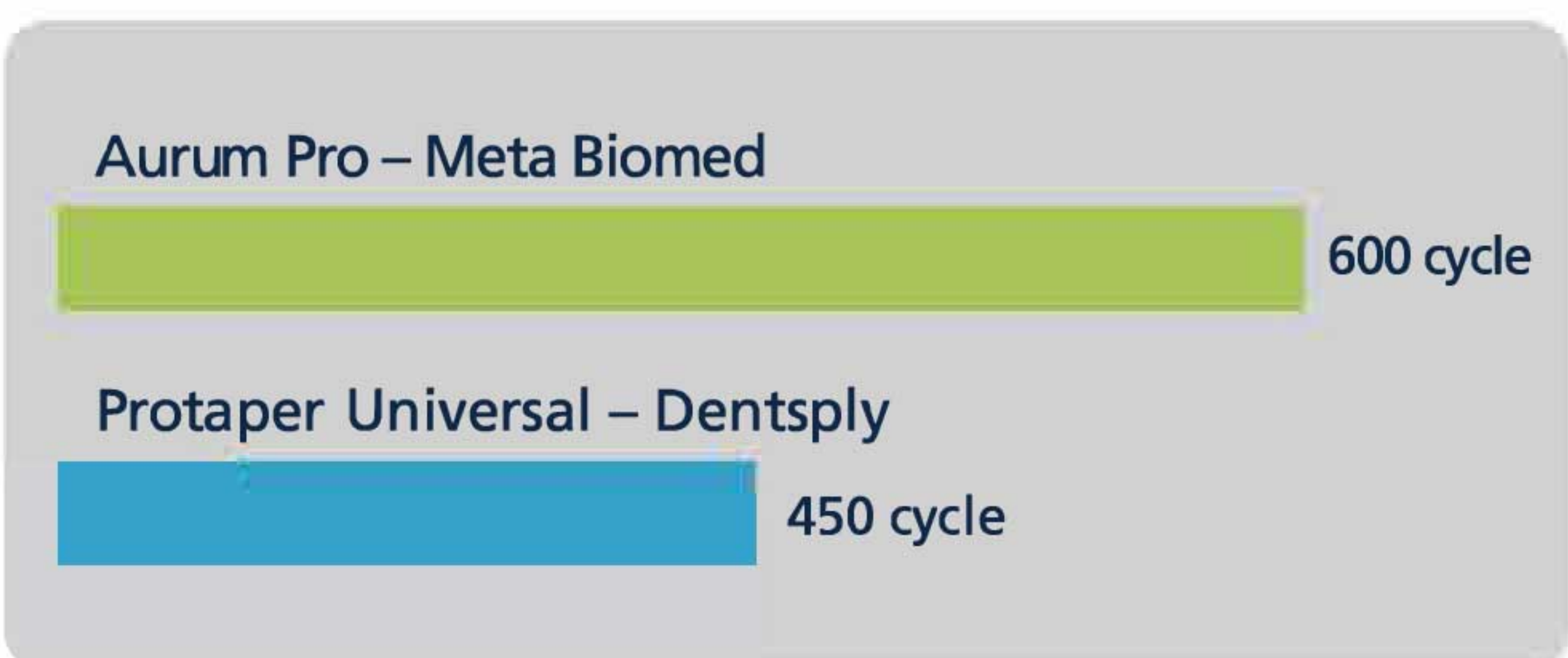Aurum Blue Triangle™

Controlled Memory File

فایل Aurum Blue Triangle™

20/10 20/04 25/06 35/04

Blue; the color of extreme elasticity & fracture resistance

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Advantage

• Martensitic Alloy

Aurum Blue files are manufactured utilizing a unique process in which the crystallographic phase transition from austenite to martensite at room temperature in contrast to conventional NiTi files, making the files extremely flexible and fracture resistant


• Super Flexible Alloy

Significantly reduce the risk of ledging, transportation or perforation of a canal

• Pre-bendable

more easier access to a canal orifice


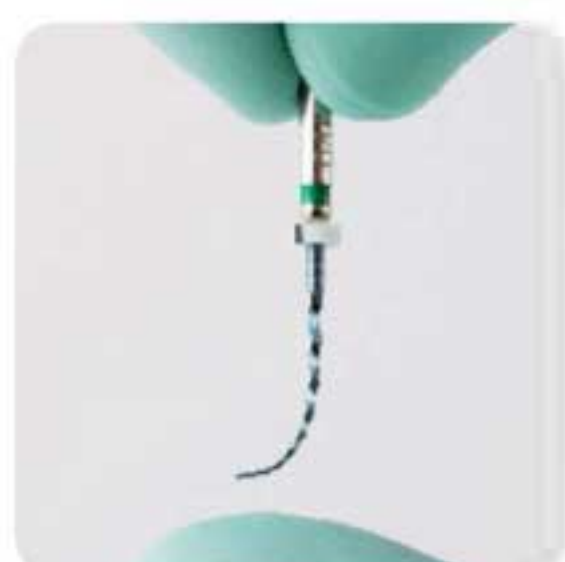
- Surface Treatment

Removing imperfection

Comparison

• Extremely resistant to fracture

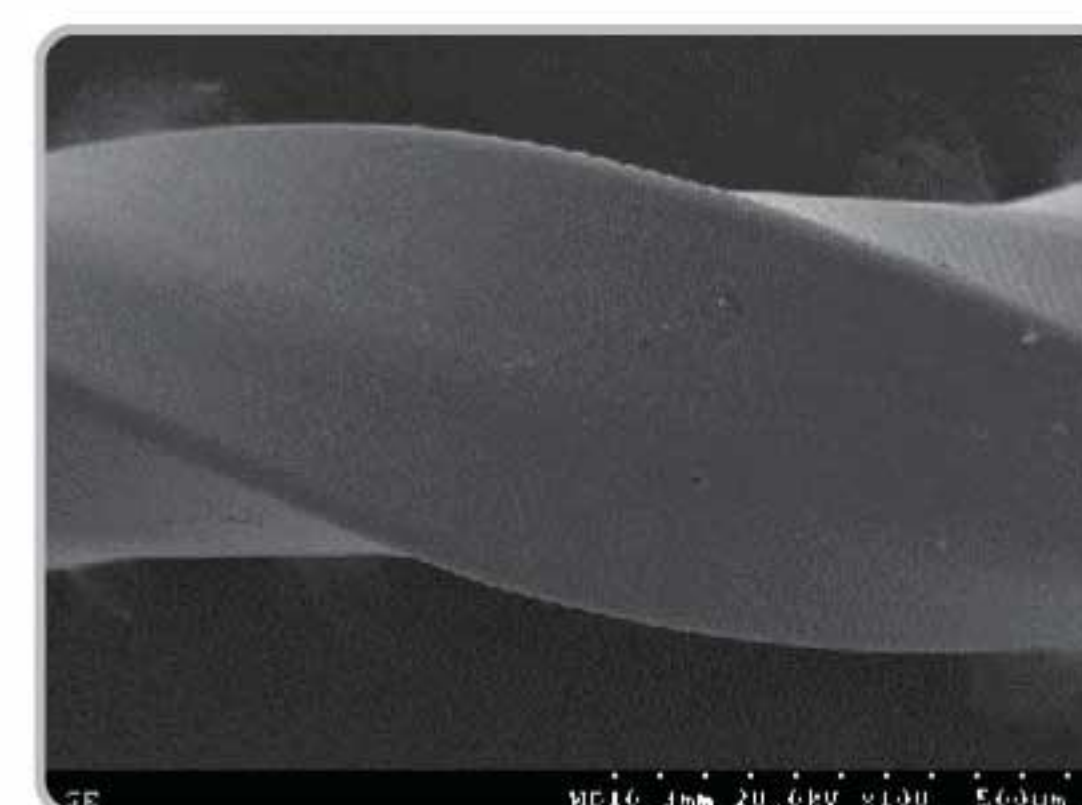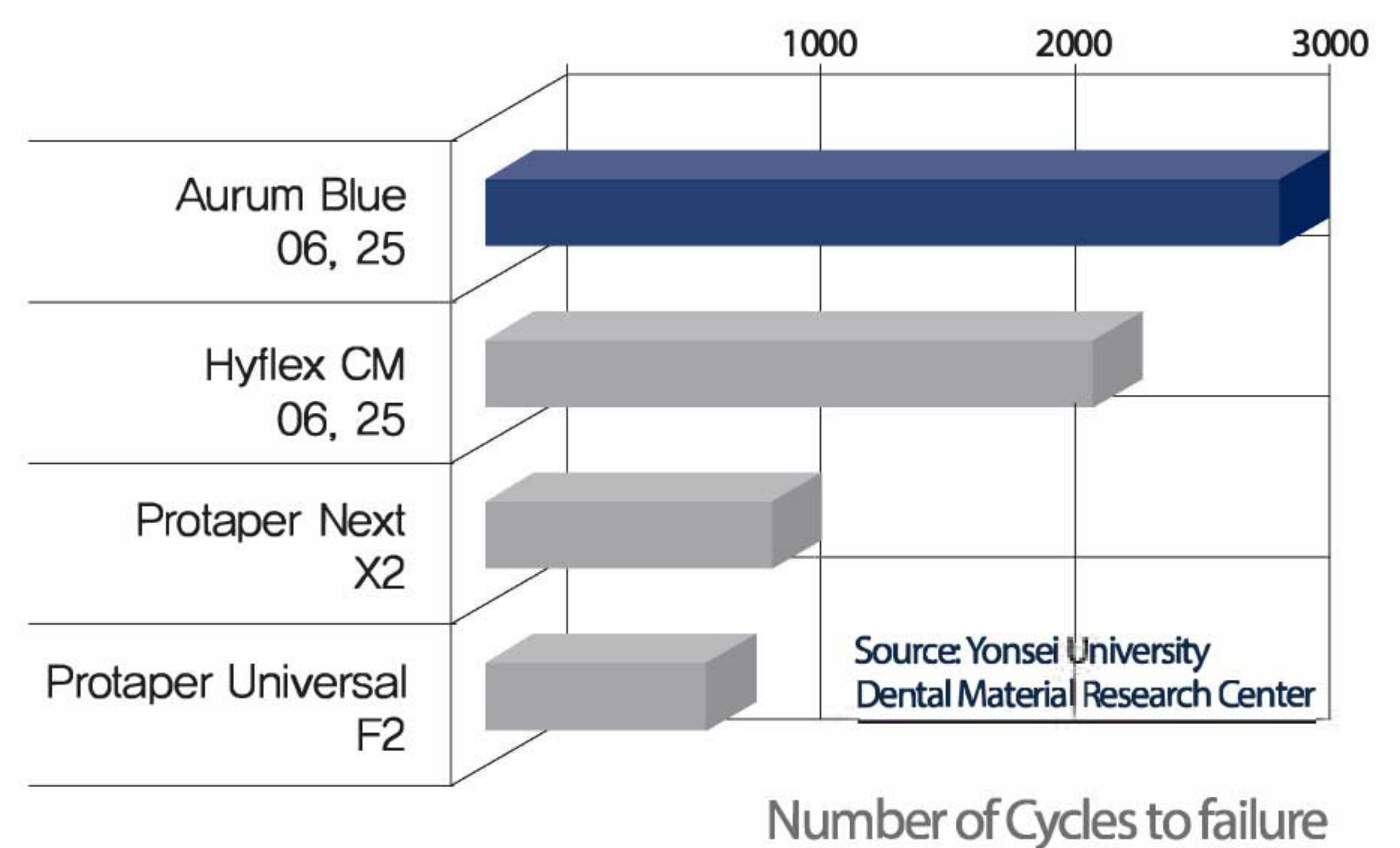• Aurum Blue

• Protaper Universal

Aurum GTM

Glide Path File

13/025 17/025

Blue; the color of extreme elasticity & fracture resistance

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
Advantage

• Increased Flexibility

With increased flexibility and more resistance to cyclic fatigue due to special heat-treatment, Aurum G offers many advantages over manual stainless steel files;

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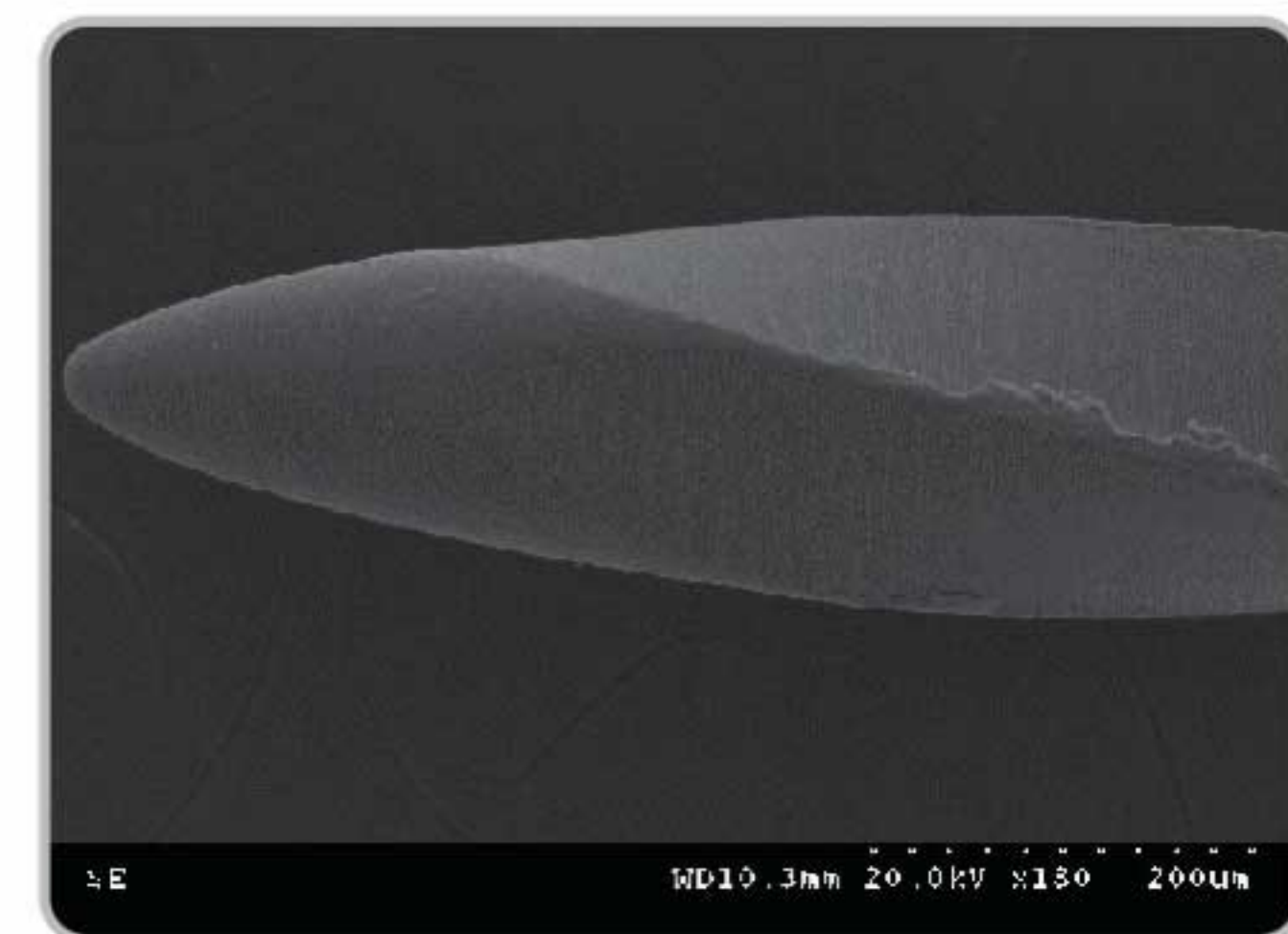
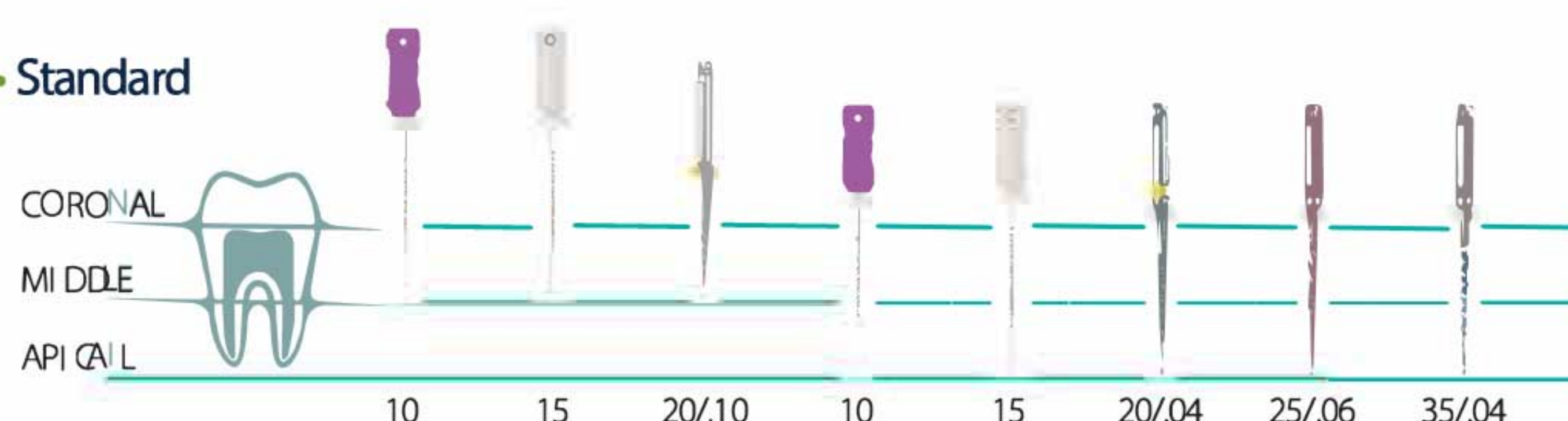
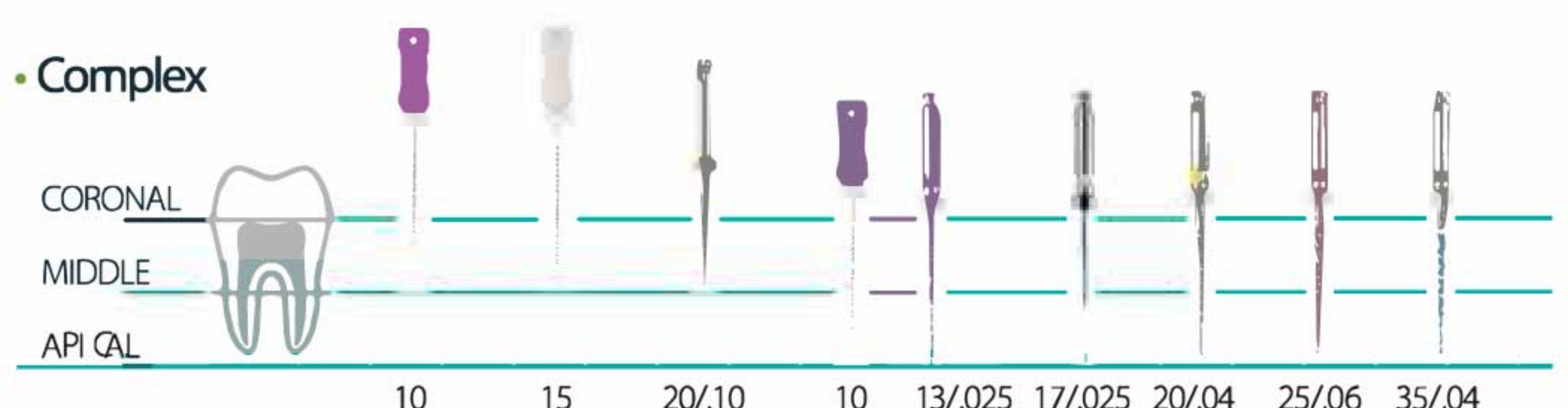
Length Marked | Patent No.: US 7097454
For Accurate Depth Measurement and Insertion

- Time saving with easy filling For (4% - 6% - 8%)
- For faster and easier root canal obturation
- Stiff but flexible to provide ideal working qualities
- Uniform and non-distorting crimp
- Long shelf life and freshness assured
- Suitable for vertical and lateral condensation
- Variety in the selection of sizes and packaging available

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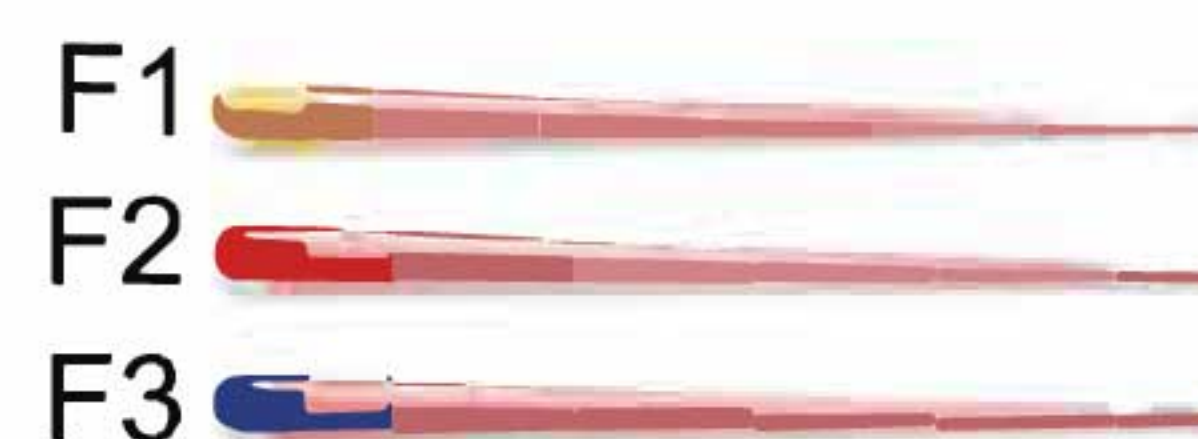Available Sizes

- ISO sizes : 10, 15, 20, 25, 30, 35, 40, 45, 50, 55, 60, 70, 80
- Non-Standardized (Accessory sizes): XF, FF, MF, F, FM, M, ML, L, XL
- Special Taper Points : .04, .06

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Advantages


- Excellent biocompatibility
- Easy to mix paste-paste
- Hermetic sealing ability
- Nonstaining to the teeth
- Insoluble in tissue fluids
- Good radiopacity

Components


- 13.5g dual syringe
Base 9g (Epoxy resin, Calcium phosphate)
Catalyst 4.5g (Amines, Bismuth subcarbonate)
- Spatula
- Mixing Plate
- Dimension : 55x180x40(mm)
- Weight : 98g

Physical and Chemical Properties


Flow

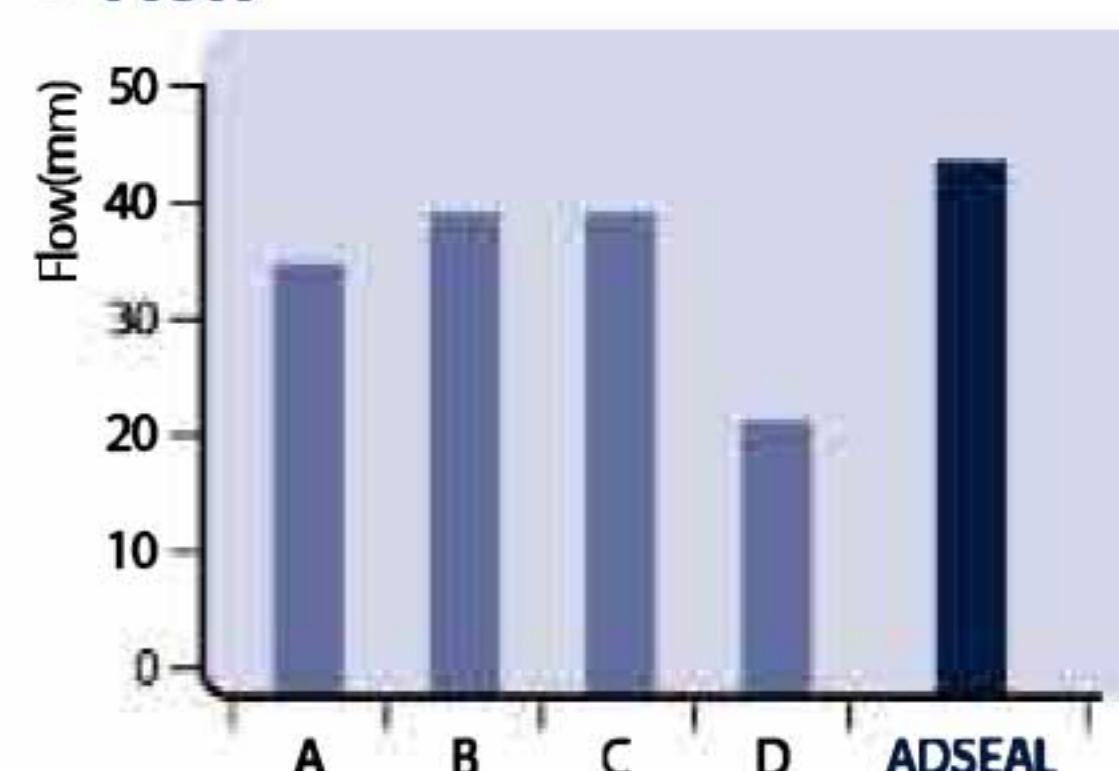
Solubility

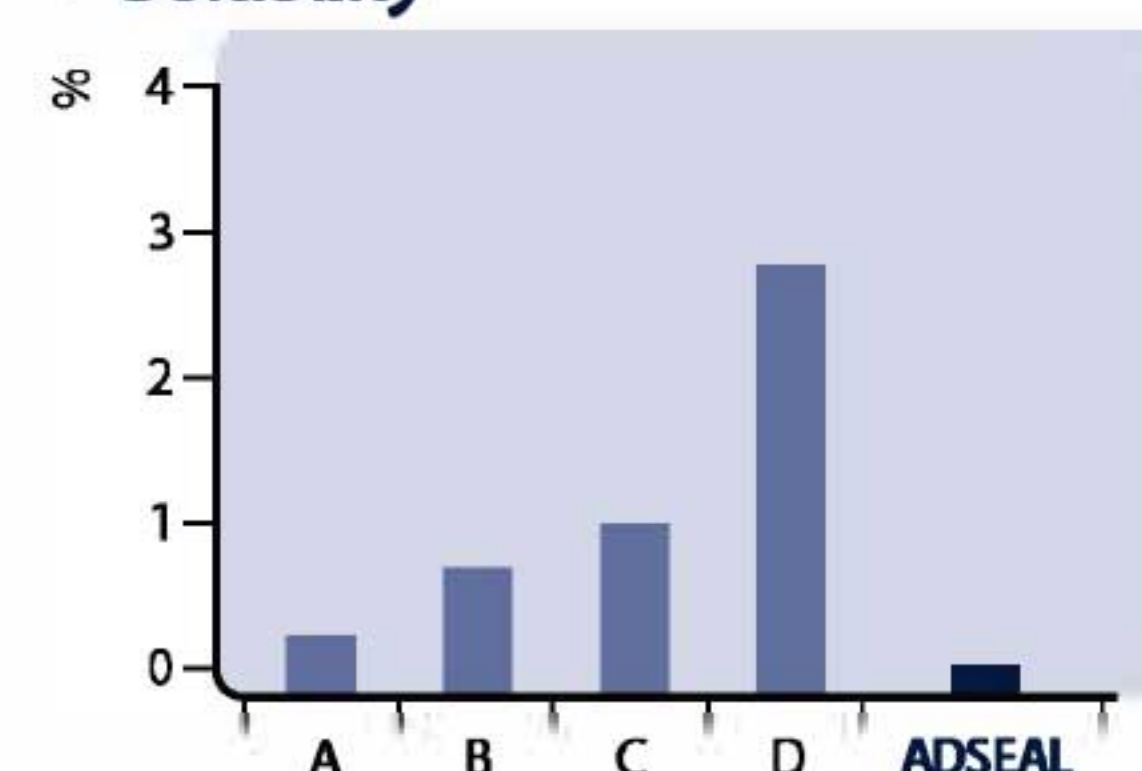Radiopacity

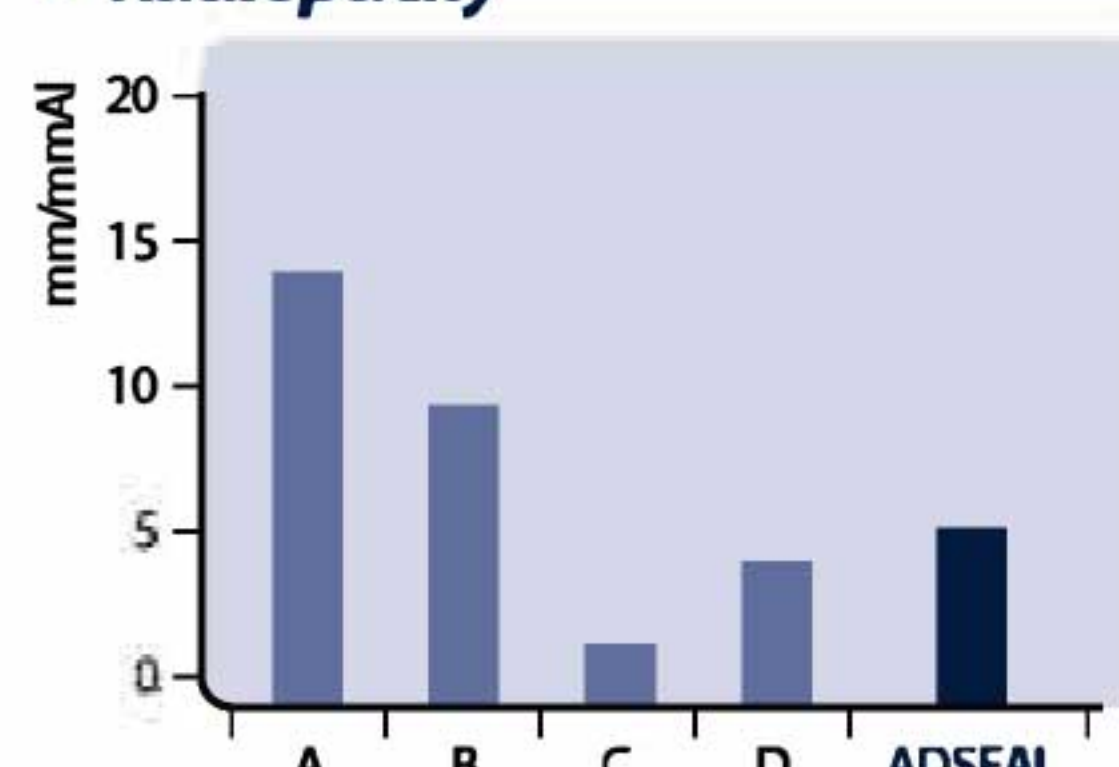Material Properties

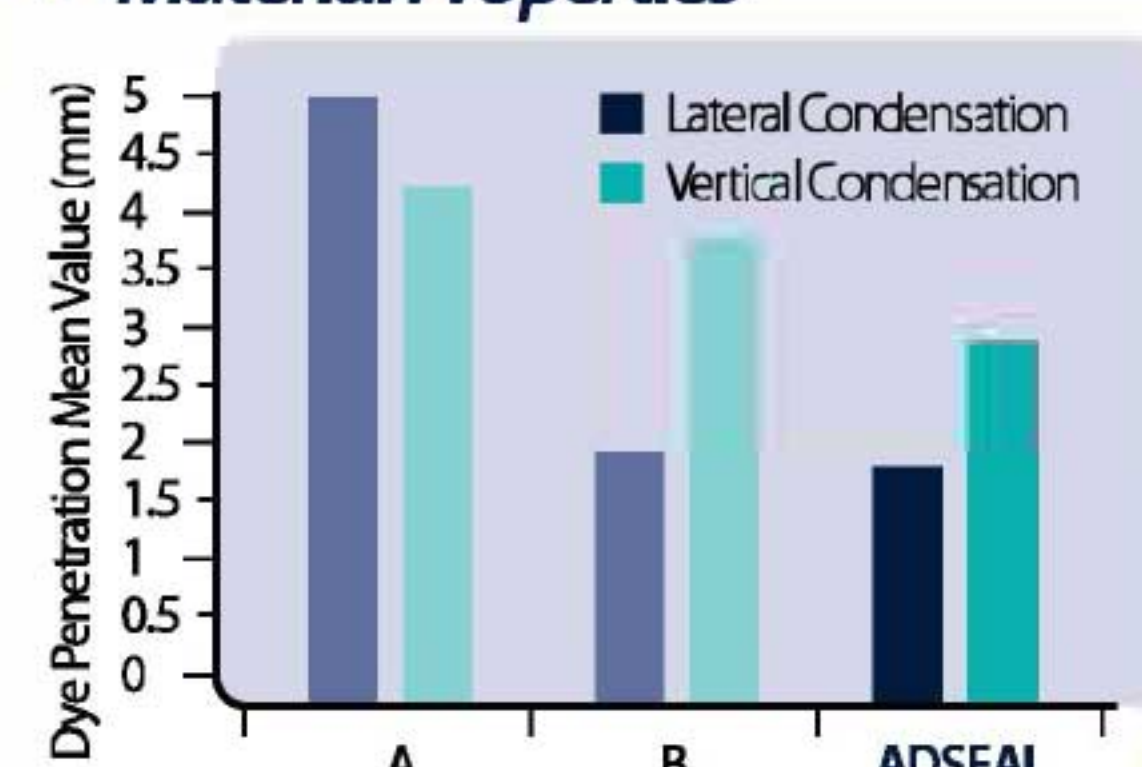- Working Time : 35 Minutes at 23°C (73°F)
- Setting Time: 45 Minutes at 37°C (99°F)
- Film Thickness: 3.3 μ m

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Applications

- Root canal preparation with files and reamers
- Presence of calcifications in the canal
- Narrow canals in posterior teeth
- Smear layer removal

Composition

- EDTA, PEG

Advantages

- Removal of inorganic substances by chelation (EDTA-Ca²⁺)
- Easier filing with lubricating action
- Produces the cleanest canal wall after smear layer removal

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Applications

- Weeping canal
- Formation of a hard tissue barrier (Apexification)
- Root canal filling for deciduous teeth
- Internal & external root resorption

Composition

- Calcium Hydroxide, Iodoform, Silicon oil

Advantages

- Excellent antibacterial effect and radiopacity
- Premixed paste
- Syringe type for easy delivery of the paste into the canal
- Disposable tips for prevention of cross-contamination

Components

- 2.2g paste in a syringe (2 syringe packing available)
- Disposable tips
- One ring rotator for direction control of the tip
- Dimension : 228x61x22(mm)
- Weight : 46g

Shelf life

- Two years from the production date.

Clinical Cases

Apexification

Pulp necrosis on #21 of 8 year-old boy took place by trauma.

About 1 year after Metapex treatment

Root Resorption

Internal resorption on disto-lingual root

Root perforation by inflammatory resorption

Metapaste

Calcium Hydroxide With Barium Sulfate

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

Physical sealing : Jamming effect

As an orthograde intra canal grafting material, OrthoMTA can achieve perfect sealing of root canal space. This is due to jamming effect of 2 micron particles at the apical foramen. OrthoMTA has superior sealing properties to any other materials.

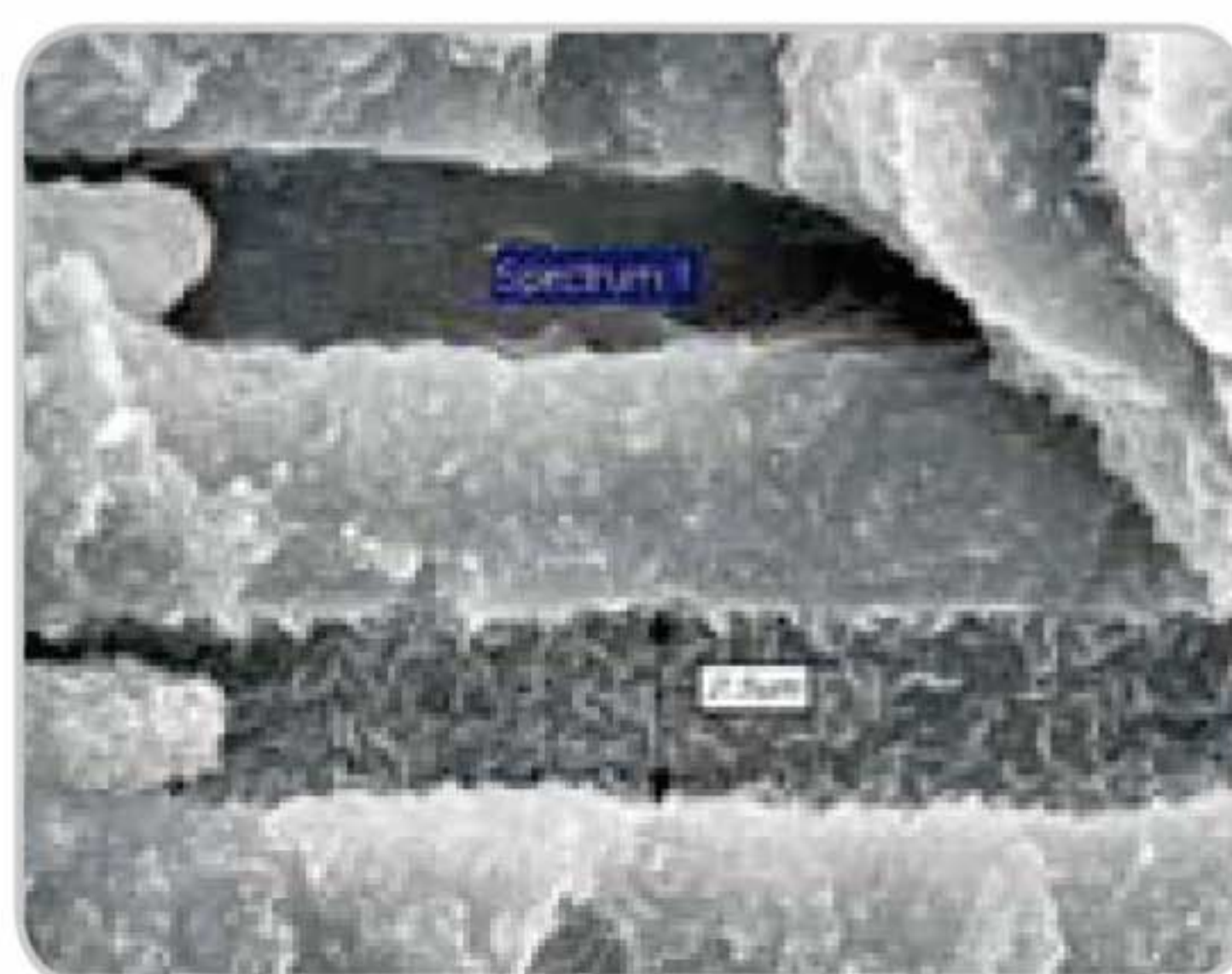Chemical sealing : C-S-H & Hap

This is due to its hydrophilic properties under moist conditions. Short tags of C-S-H obturate the opening of the dentinal tubules. As a result, the entombing effect can be induced. Long tags of hydroxyapatite fill the space of it and achieve fossilization of bacterium.

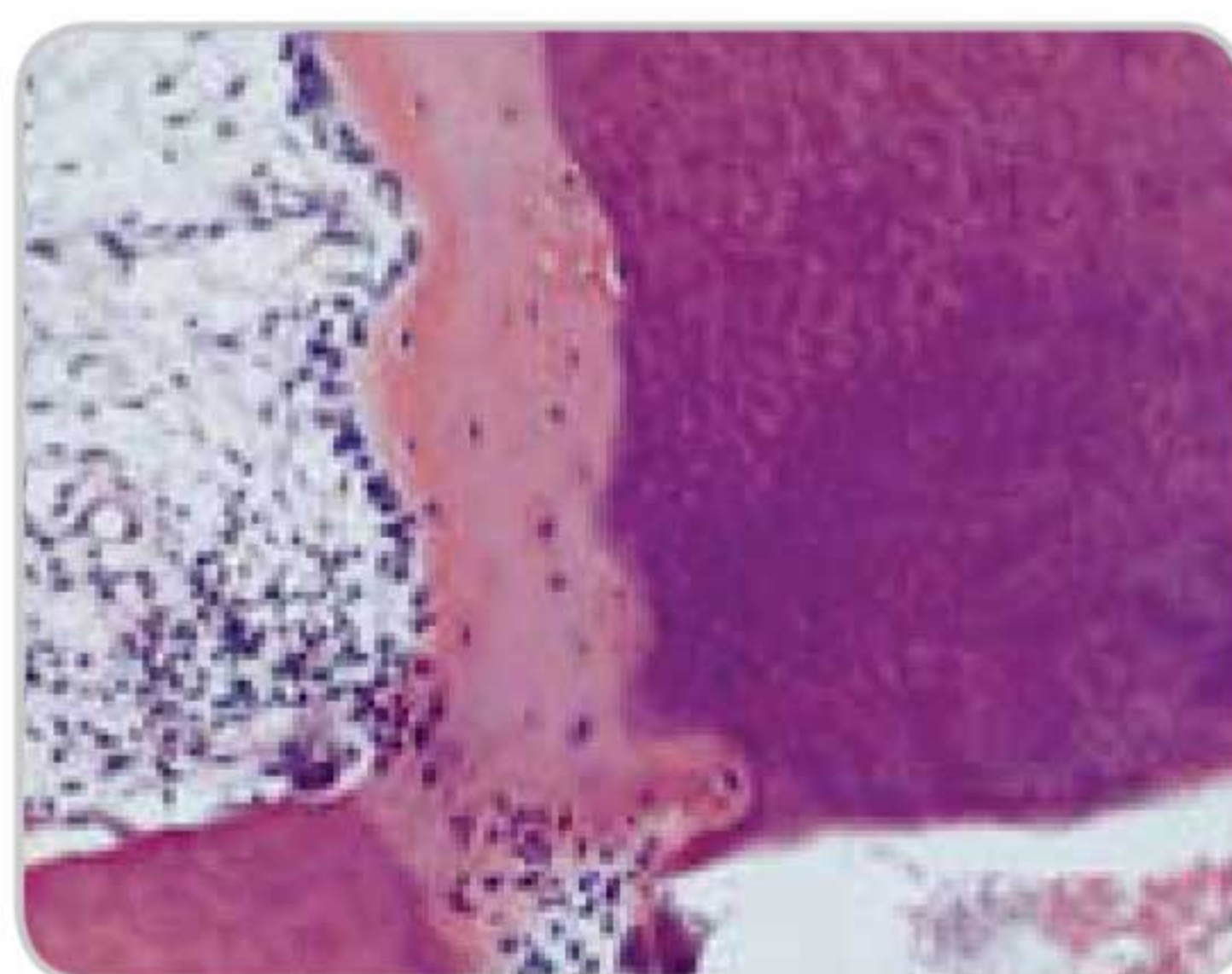Biological sealing : New cementum

New cementum is regenerated at the apex of root after orthograde intra canal graft. Because OrthoMTA is bioactive material, damaged periodontium is regenerated. Biological sealing can be induced by OrthoMTA.

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
1. 0.3g RetroMTA powder and wetting solution

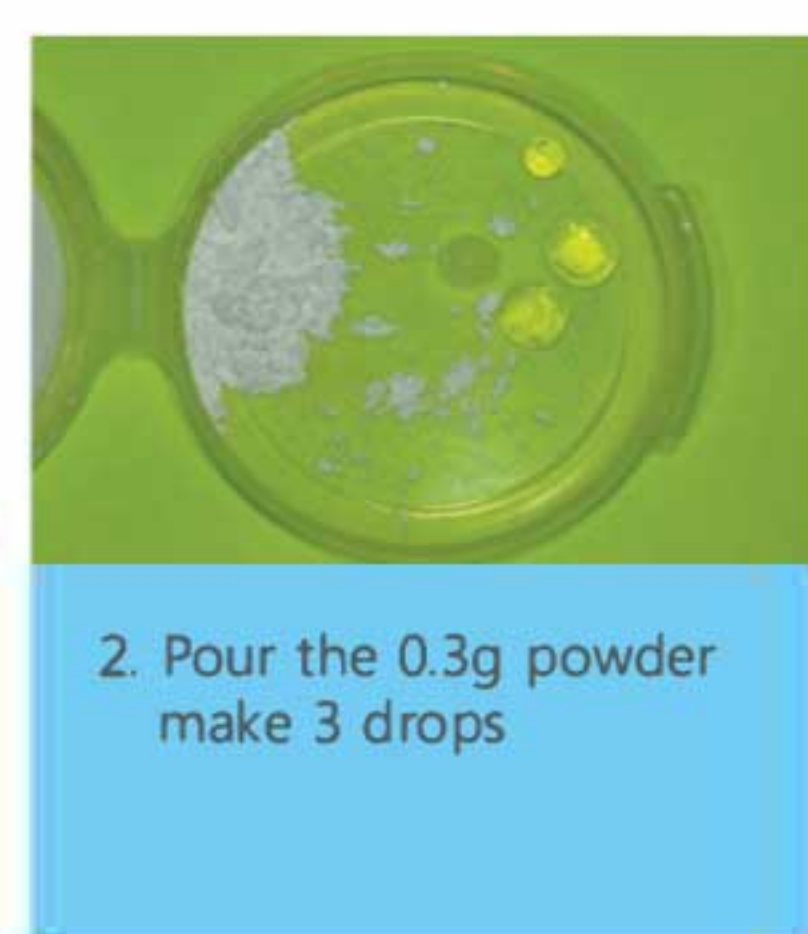2. Pour the 0.3g powder make 3 drops

3. Wet it gently within 20 seconds

4. Wait until the shiny surface disappears

5. Apply it with a carrier

Composition and Ingredient Information

Product	Chemical Name	CAS No.	Content (wt%)
RetroMTA (OrthoMTA II)	Calcium Carbonate (CaCO_3)	471-34-1	60~80
	Silicon Dioxide (SiO_2)	7631-86-9	5~15
	Aluminum oxide	1344-28-1	5~10
	Calcium zirconia complex	contrast media	20~30
	Total		100

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DX-Condenser

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Technical Specifications

- DX Condenser 35 / 70 (Green, 0.35mm Ni-Ti tip & 0.7mm Stainless Steel tip)
- DX Condenser 40 / 80 (Black, 0.4mm Ni-Ti tip & 0.8mm Stainless Steel tip)
- DX Condenser 50 / 100 (Yellow, 0.5mm Ni-Ti tip & 1.0mm Stainless Steel tip)

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EsCom 100

Light-Cured Restorative Composite Resin

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Applications

- Class I , II ,III,IV, V restorations
- Diastema closures
- Direct veneers
- Composite and porcelain repairs

Composition

- Bis-Gma, UDMA, Bis-Ema
- Barium aluminium boro silicate


Advantages

- Low 1.6% shrinkage for superior marginal integrity
- Various shade
- Filled with 40 nano sized filler for excellent polishability
- Optimum handling - Excellent Sculptability for restoring natural shape and function.
- Non-Slumping. Non-Sticky

Components

- Nexcomp Kit
 - 4g x 7 Syringe (A1, A2, A3, A3.5, B2, A2E, B2E)
 - Meta P&Bond 5g
 - Meta Etchant 3g
- Nexcomp Refill
 - 4g x 1 Syringe (A2E, B2E, A1, A2, A3, A3.5, A4, B1, B2, B3, C2, C3, D2, D3, TL, WT, OP, OA2, OA3)g

Physical and Chemical Properties

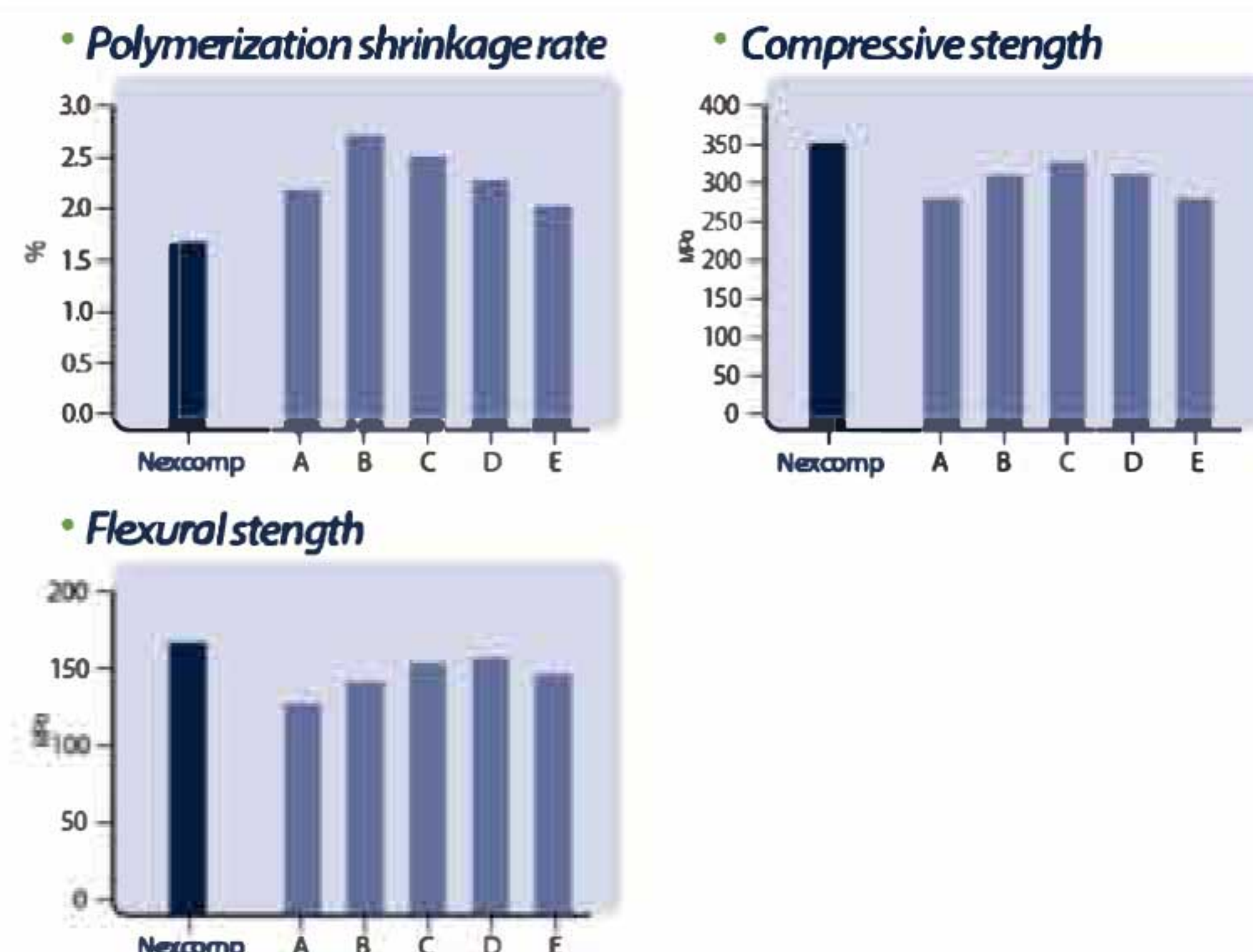Clinical Cases

Shelf life

- Three years from the production date.

EsFlow

Light-Curing Flowable Composite Resin

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Highly aesthetic, nanohybrid flowable composite.
Medium viscosity allows gentle flow for easy handling and precision application.

- Contains high molecular weighted urethane methacrylate monomer for superior compressive strength
- Light-cured to a hardness equivalent to natural enamel
- Highly radiopaque

Syringe pkg

Contains : 2 x 2g Syringe, 8 Dispensing Tips

Specify :

#3121A1	(A1)	#3121B1	(B1)
#3121A2	(A2)	#3121B2	(B2)
#3121A3	(A3)	#3121BL	(Blue)
#3121A35	(A3.5)		

Bent Dispensing Tips

Pkg of 5,000 or 10,000

#911120 20 ga. (black)

Contains : 20 ga. dispensing tips

Nexcomp Flow

Nano Hybrid Flowable Composite Resin

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Applications

- Restoration of invasive cavity
- Pits and fissure sealing
- Base and Liner for direct restoration
- Repair acrylic temporary materials
- Repair of esthetic indirect restorations

Composition

- Bis-Gma, UDMA, Bis-Ema TMPTMA
- Barium aluminium boro silicate

Advantages


- Low viscosity, light-cured, radiopaque flowable composite resin for minimally invasive cavity preparation.
 - Filled with 40nano sized filler for excellent polishability
 - Flow on demand handling
 - Excellent Radiopacity
 - Available in 14 shades including opaque and translucent
 - No slumping and running

Components


- Nexcomp Flow set
 - 2g X 5 syringes(A2, A3, A3.5, B2, B3), 20ea disposable tip
- Nexcomp Flow set
 - 2g X 2 syringes(A1, A2, A3, A3.5, A4,B1, B2, B3, OA2, OA3, C2, D2, WT, TL), 10ea disposable tip

Physical and Chemical Properties

• Compressive Strength(MPa)

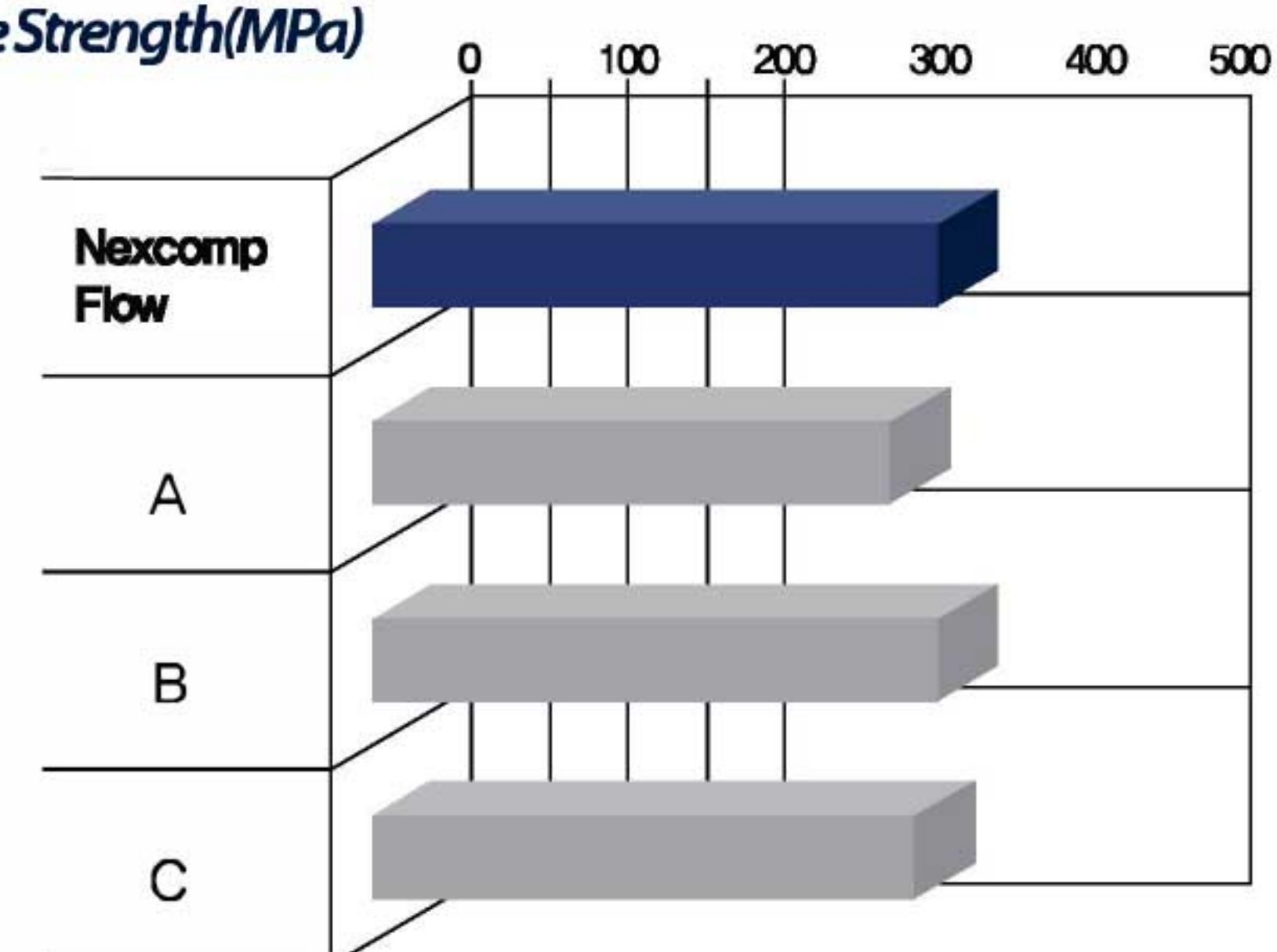• Flexural Strength(MPa)

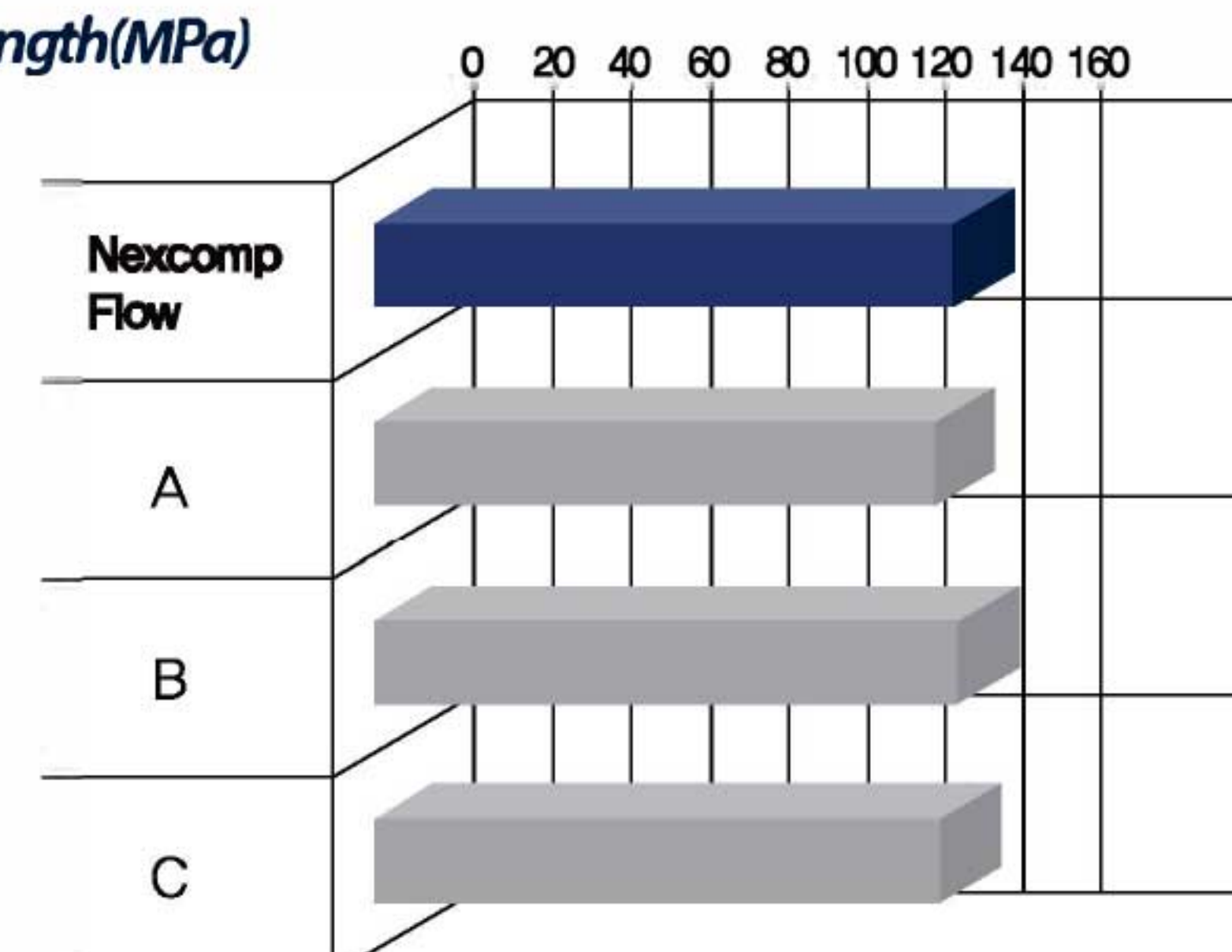Shelf life

- Three years from the production date.

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Contents

Core .it Kit : Core .it Dual 10g x 2syringes
FineEtch 5ml x 1syringe, 4Tips, EsBond 5ml x 1bottle,
Activator 5ml x 1bottle,Tip Holder, Dispensing Dish, Micro Tip Application 30ea

Contents

Cartridge type (Automix) : 1ea(Base 25g + Catalyst 25g)

NexCore

Dual-Cured core build-up composite Resin

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Applications

- Post adhesions and core build-up in restorations

Advantages

- Suitable viscosity and flow for post cementation and core build-up
- Dual curing: Light-curing and self-curing are possible so there is no need to worry about the polymerization failure of resin
- Biocompatibility (ISO10993)
- Easy to mix by auto-mixing type, Direct injection: Intraoral-Tip, Endo-Tip
- Radiopaque
- 3 Shade (Tooth, White: Ceramic restoration Blue: Metal restoration)
- Cuttability likes dentin (Flexural modulus: NexCore 14.2GPa)
- Three years from the production date

Components

- NexCore Cartridge 50g paste 1ea, Automix-tip 30ea, Intraoral-tip 30ea
- NexCore Cartridge Refill 50g paste 1ea (Shade: Tooth, White, Blue)
- NexCore Dual Syringe 9g paste 2ea, Automix-tip 20ea, Intraoral-tip 10ea, Endo-tip 10ea
- NexCore Dual Syringe Refill 9g paste 1ea (Shade: Tooth, White, Blue)


Physical Properties

- Composition
Barium glass, Ytterbium trifluoride, fumed silica in a Bis-GMA based matrix of dental resins
- Physical Properties


Working time	1:50 ~ 2:20 min
Setting time	4:00 ~ 5:00 min (Self-cure) 20 ~ 40 sec (Light-cure)
Flexural strength (light curing)	115.52 MPa
Water Sorption, Solubility	16.5/ 2.3 $\mu\text{g}/\text{mm}^2$
Radio-opacity	Optical density was equivalent with optical density at 4.2mm aluminum step wedge
Shade and color stability	No change in color and identical to the shade guide.

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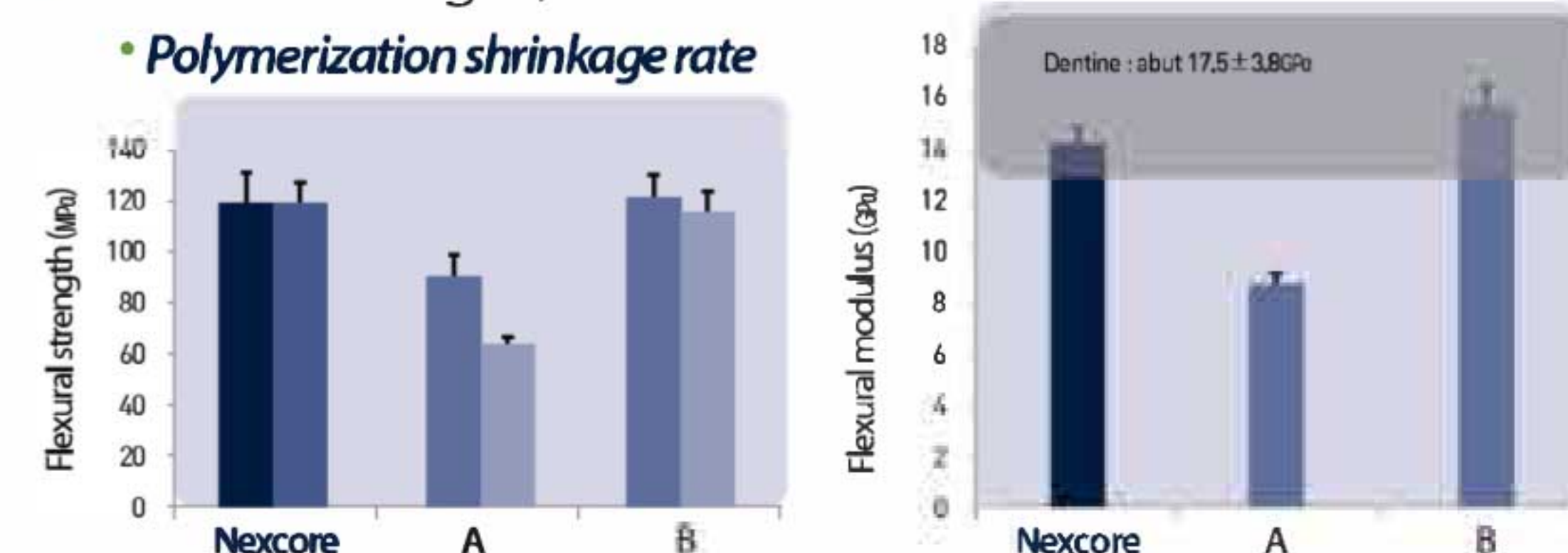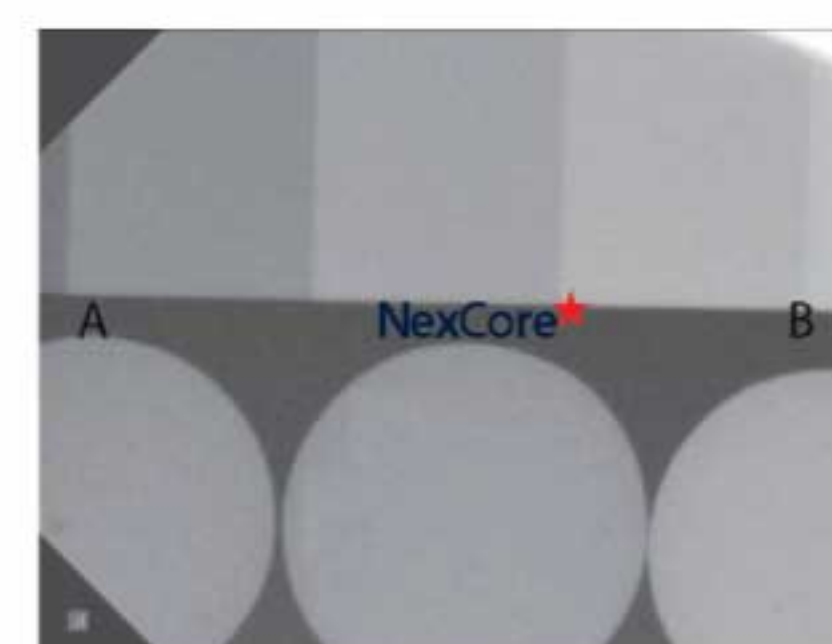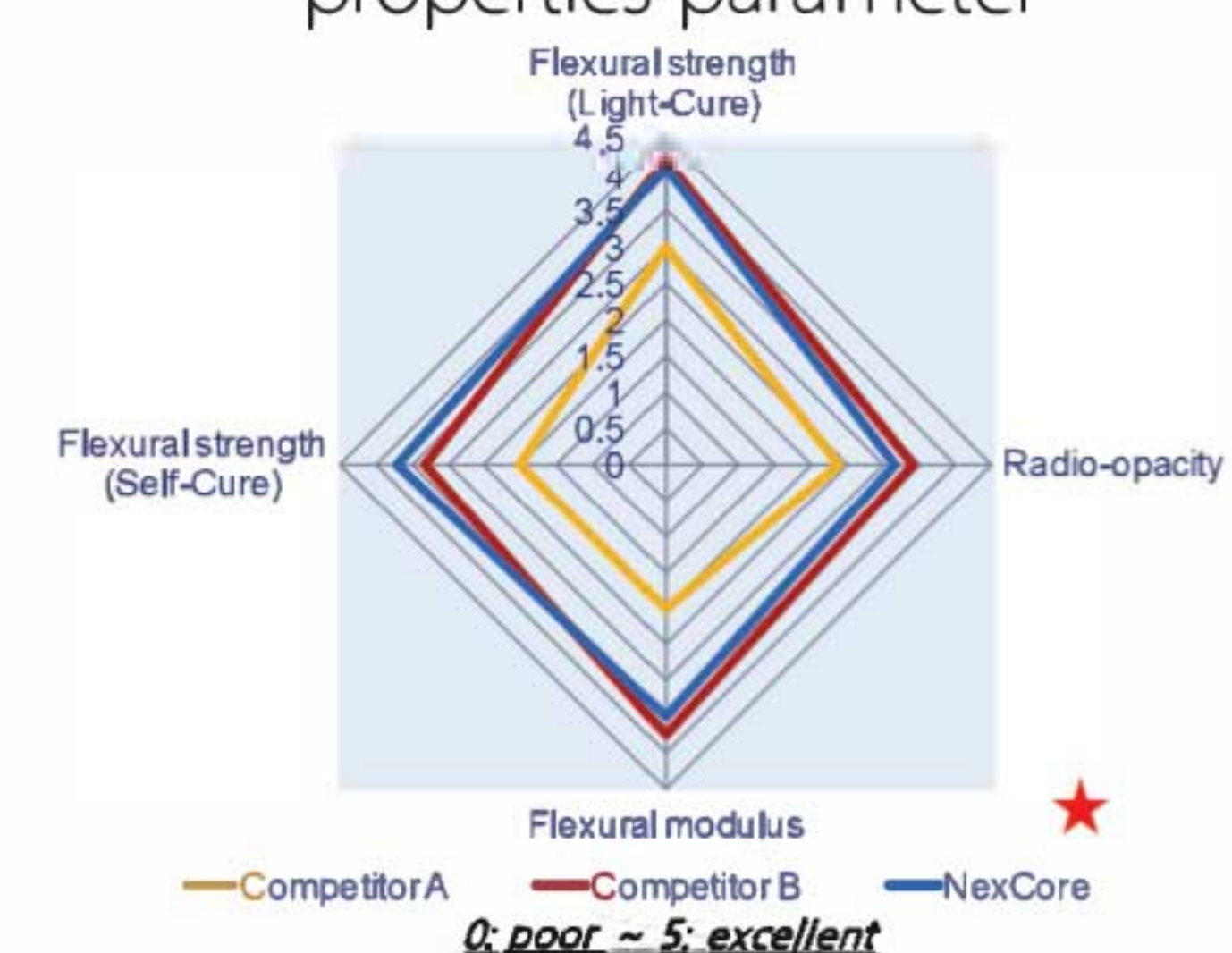
Applications

- Crowns & Bridges (ceramic, composite, porcelain)
- Inlay/Onlay cementation
- Bonding for porcelain veneers

Composition

- Base-Bis-Gma, TEGDMA,
Barium aluminium boro silicate
- Catalyst-Bis-Gma, TEGDMA,
Barium aluminium boro silicate

Advantages

- Easy to use for cosmetic indirect restorations.
- Provides high resistance to compressive force.
- Provides a dual-cure for highly shaded, thick or opaque cementation.
- 5 shades correspond with shade guide.
[A1, A3, B2, TL(Translucent white), OP(Opaque)]

Components

1. Metacem set


- 2 Dual syringe(A3, Translucent white shade, 9g/ea)
- P&Bond 5g
- 20 Mixing tips
- Spatula
- 10 Adhesive brush
- Etchant 3g
- 20 Etchant disposable tips
- Mixing pad paper

2. Metacem Refill

- 1 Dual syringe(A1, A3, B2, Translucent white, Opaque, 9g/ea)


Physical and Chemical Properties

• Flexural Strength

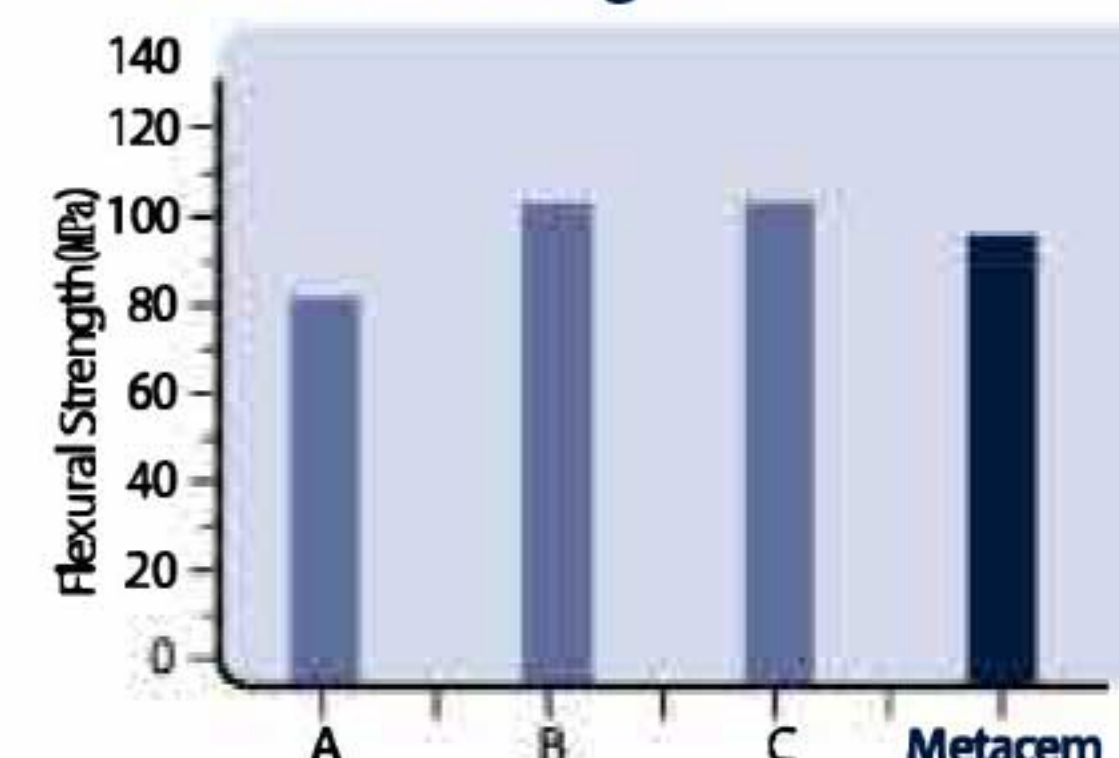
Working time 2 min
Setting time 5 min
Flexural strength 99 MPa
Radiopacity 2.33mm Al/mm

• Radiopacity

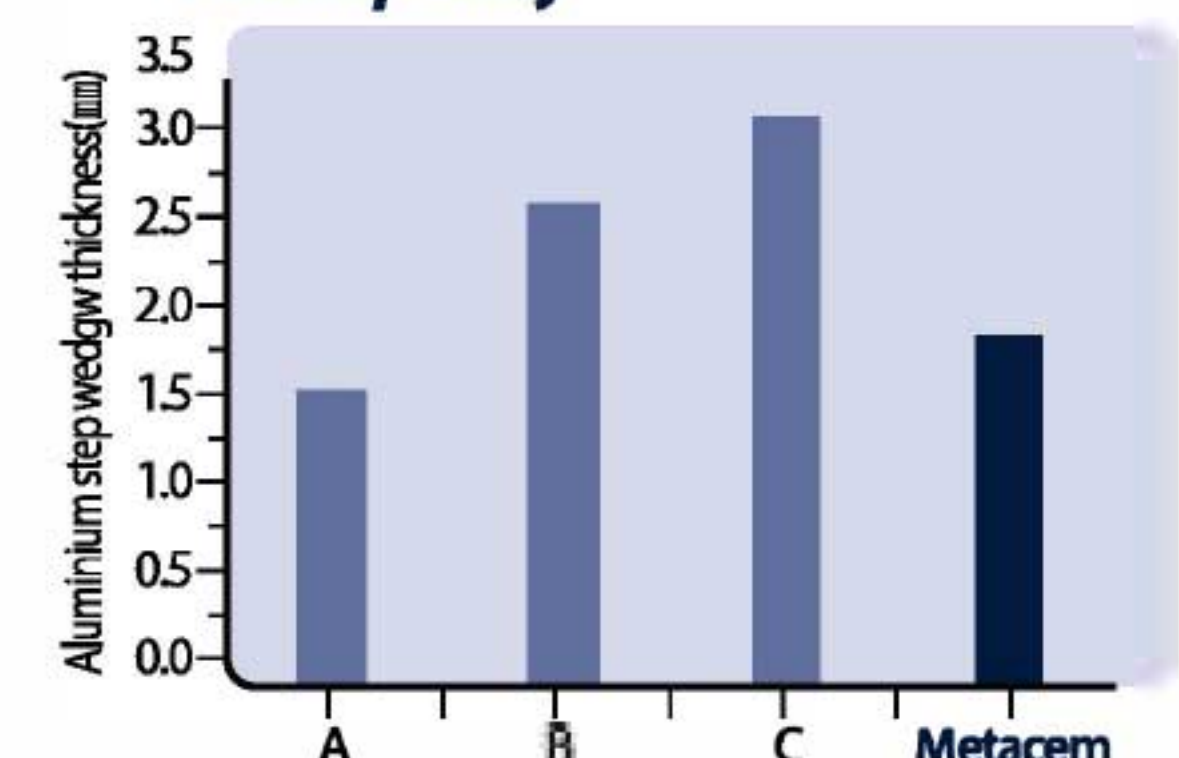Film thickness <21 μm
Water absorption 29.3 μg/mm²
Water solubility 2.3 μg/mm²

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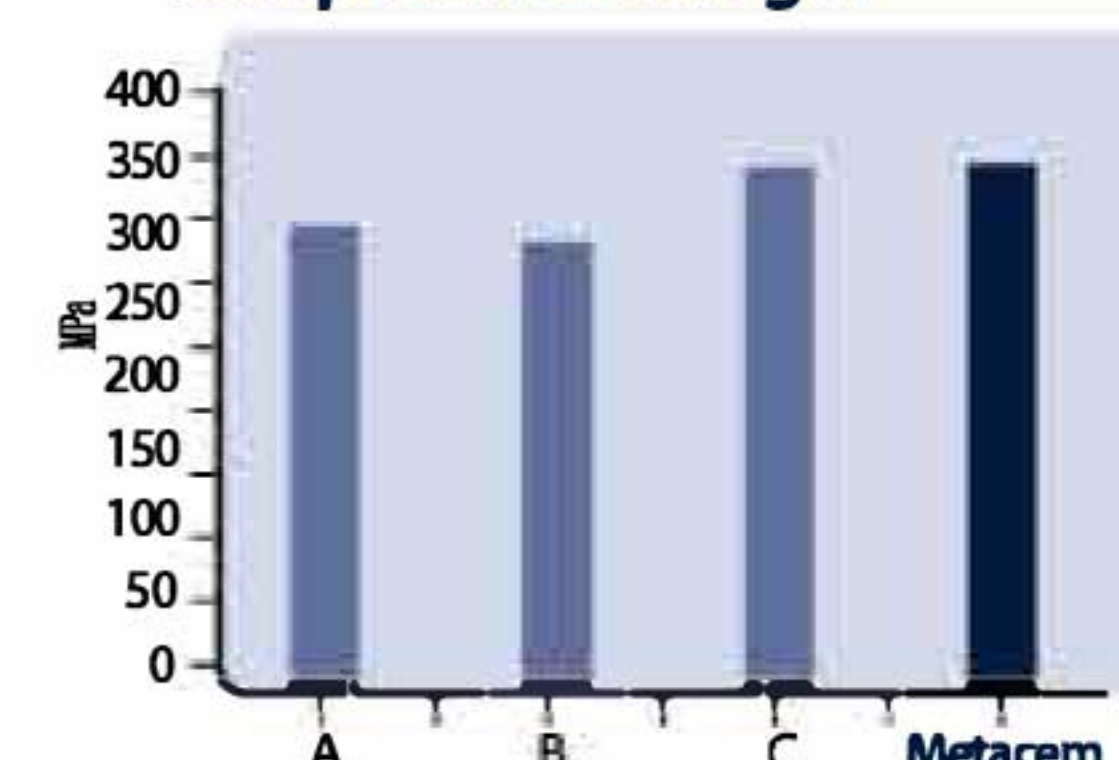Applications

- Bonding for direct composite
- Bonding for porcelain, composite, metal repair, set amalgam, post
- Bonding for porcelain veneer
- Bonding for amalgam and indirect bonding procedures with Metacem adhesive resin cement

Composition

- Bis-Gma, PMGDM, 2-Hydroxyethylmethacrylate, Ethyl alcohol

Advantages

- Superior bonding strength
- 5th Generations Adhesive
- One step priming and bonding
- Ethanol base adhesive
- Wet bonding type

Physical Properties

- Dentin shear bonding strength 10.75MPa
- Enamel shear bonding strength 30.76MPa

Components

- Meta P&Bond 5g
- Micro brush 10EA

Shelf life

- Two years from the production date.

EsBond

Light-Curing Temporary Filling Material

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
EsBond is 5th generation Bonding agent that provides high bonding strengths in a convenient single bottle application. That single component bonding agent is designed to bond all classes of direct composite restorations to dentin, enamel, as well as composite, treated porcelain and set amalgam.

1. Apply EsBond

2. Air dry

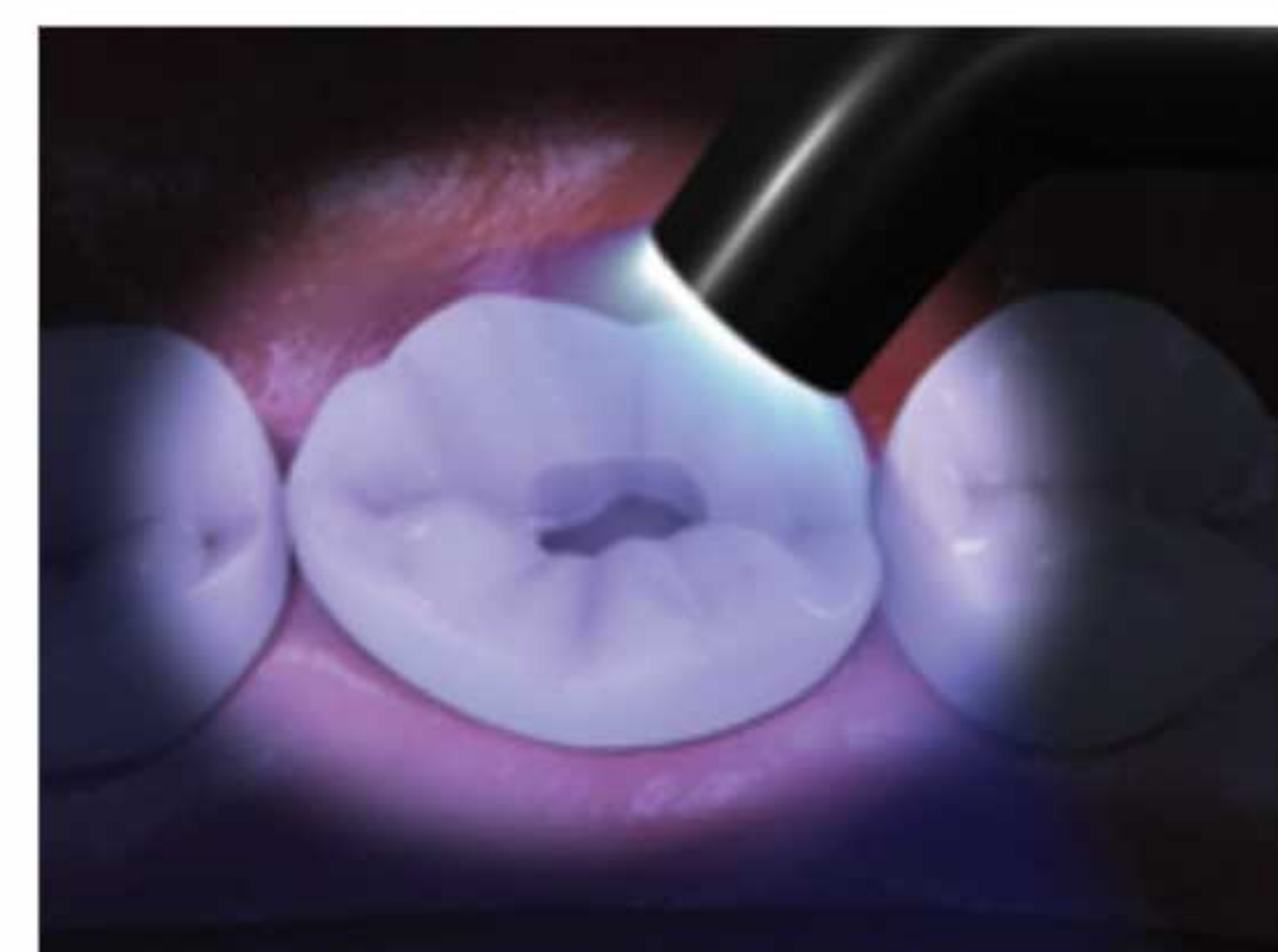
3. Light-Cure for 20 seconds

4. Place the EsFlow or EsCom100

Physical Properties

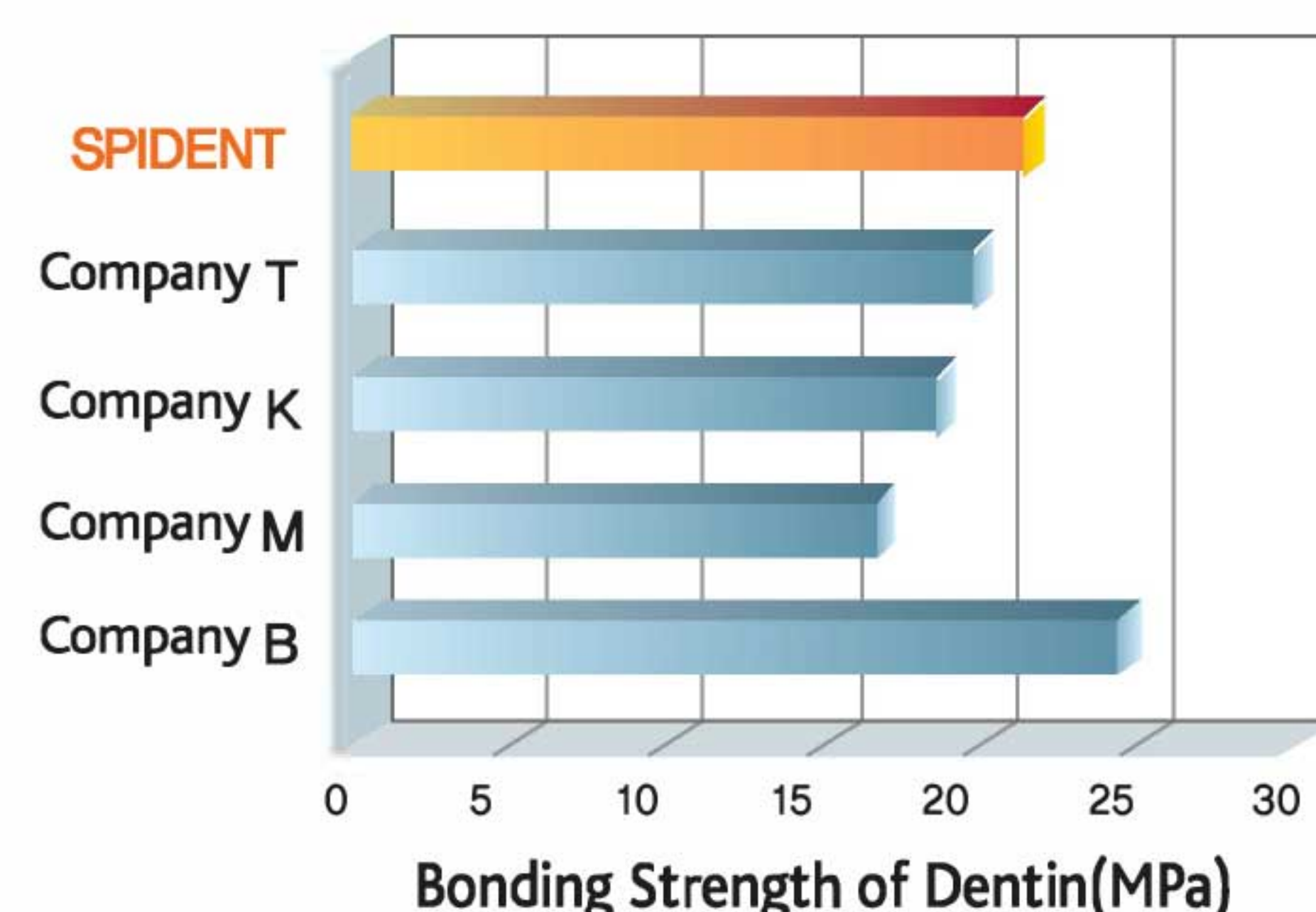Meta Etchant

Meta Etchant contains 37% phosphoric acid water thickened with unique polymer

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Seal . it is a light-curing dental material designed for sealing the enamel pit & fissure of teeth. It has the excellent characteristics of the bonding strength, low polymerization shrinkage and abrasion resistance, suitable viscosity with the Nano filler.

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EsTemp Implant is a kind of temporary resin cements and is designed to be suitable for various applications such as cementing for temporary crowns, bridges, inlays, onlays and splints.

Application

1. Implant operation

2. Apply EsTemp Implant and place temporary crown


3. Check occlusion and then remove excess material

4. Finish the restoration

Physical Properties

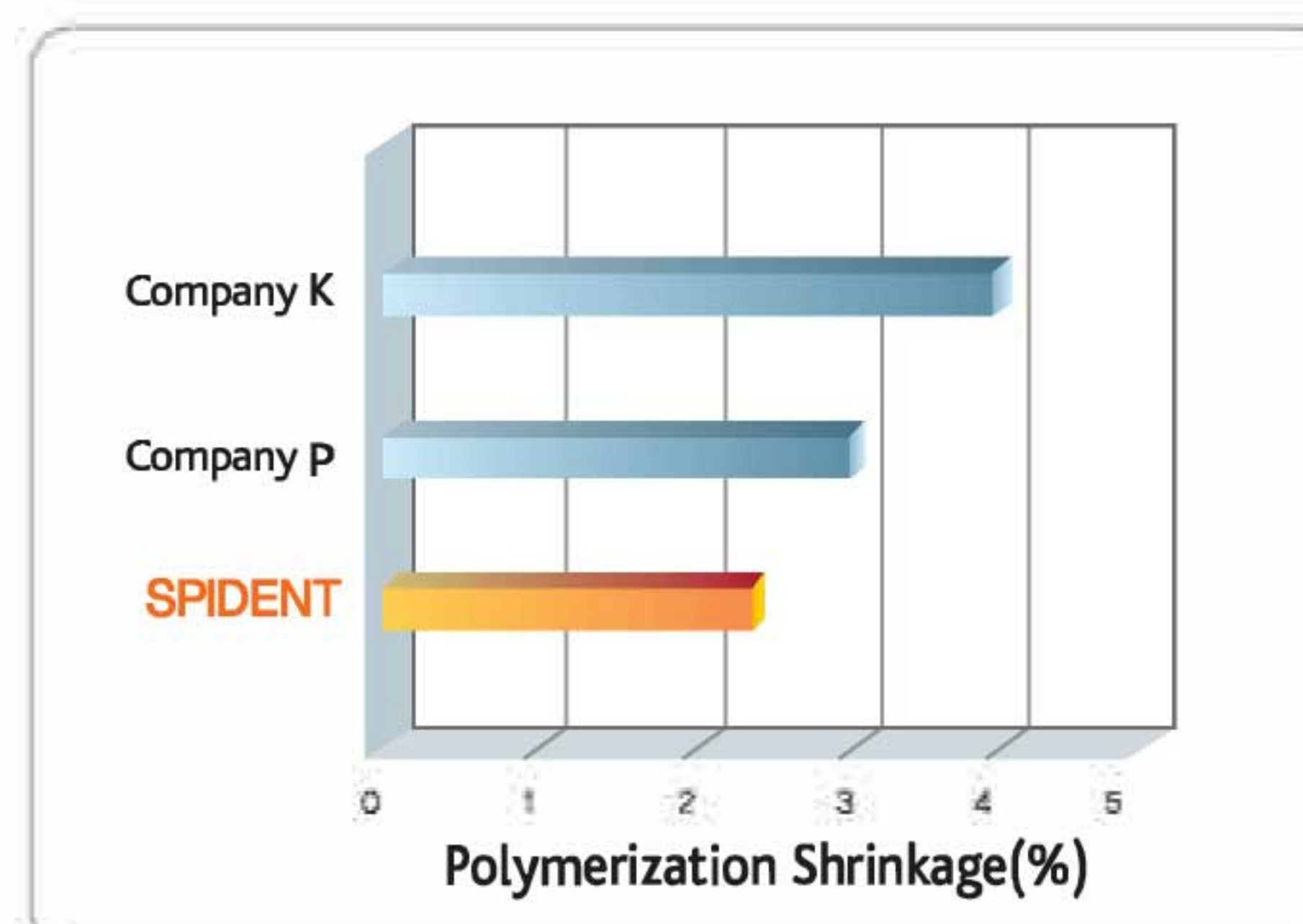Application of Auto mixing Tip

Temp.it

Light-Curing Temporary Filling Material

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NexTemp HV / NexTemp LC

A Light-Cured temporary Filling materials

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Applications

- Temporary restoration of prepared cavity(1~6days)
- Temporary filling sealing of cavities

Advantages

- High elasticity
- Easy application
- Tightness for margin sealing
- Easy removal

Component

- NexTemp HV : 3.0g Syringe x 3ea
- NexTemp LC : Syringe 2.0g paste 3ea, Disposable Tip 10ea

Shade

- NexTemp HV : TL, Blue
- NexTemp LC : TL, Blue, Yellow

Shelf life

- Nextemp HV : Two years from the production date
- Nextemp LC : Three years from the production date

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
Advantages

- Radiopaque
- Translucent
- Optimal flexural strength(close to dentin)
- Excellent light translucency for efficient light cure deep into the root canal


Composition

- Glass Fiber
- Bis-GMA based composite resin

Strength test

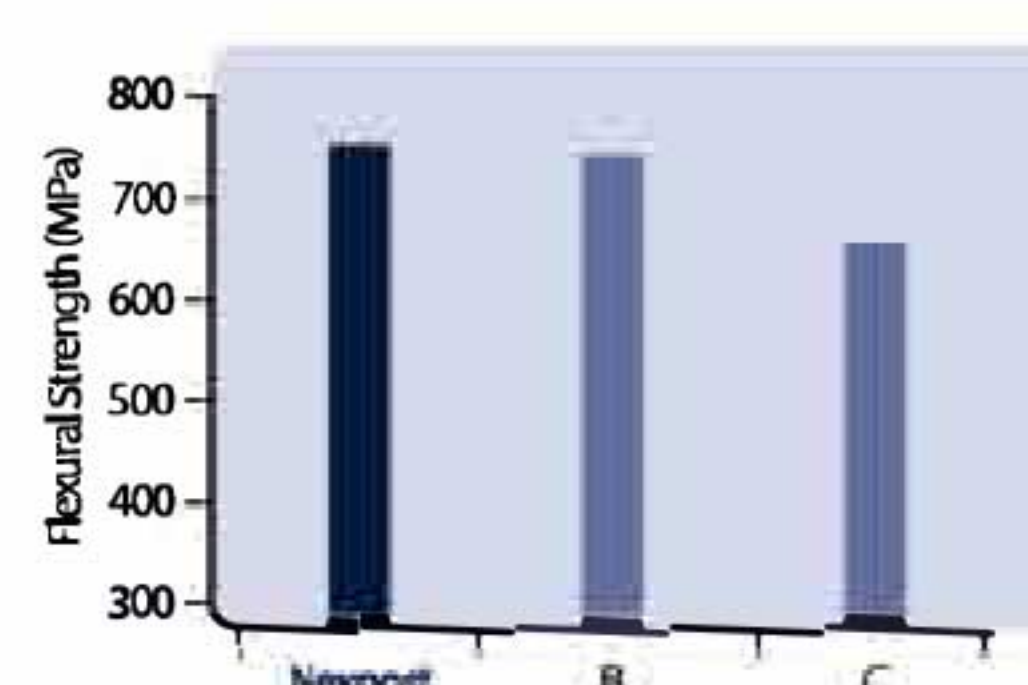
Sizes

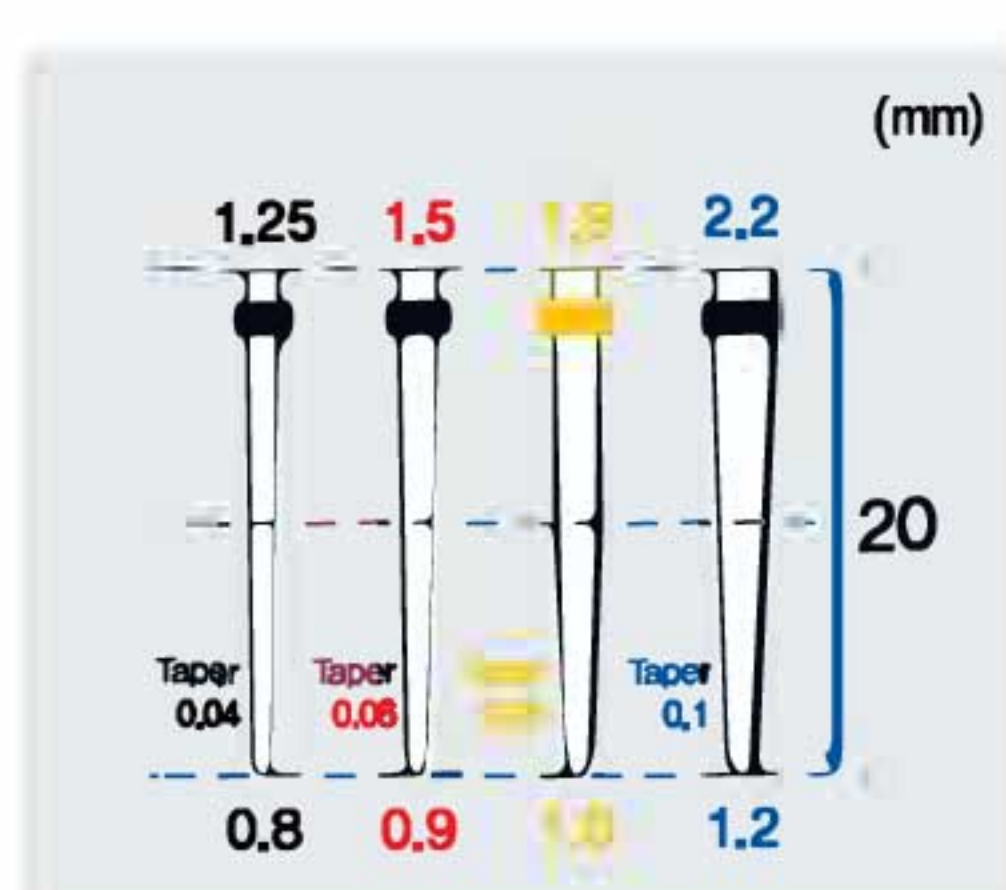
Components

- NexPost Kit (Post x 20, Drill x 5)
- NexPost Refill (Post x 5)
- NexPost Drill Refill (Drill x 1)

Physical and Chemical Properties



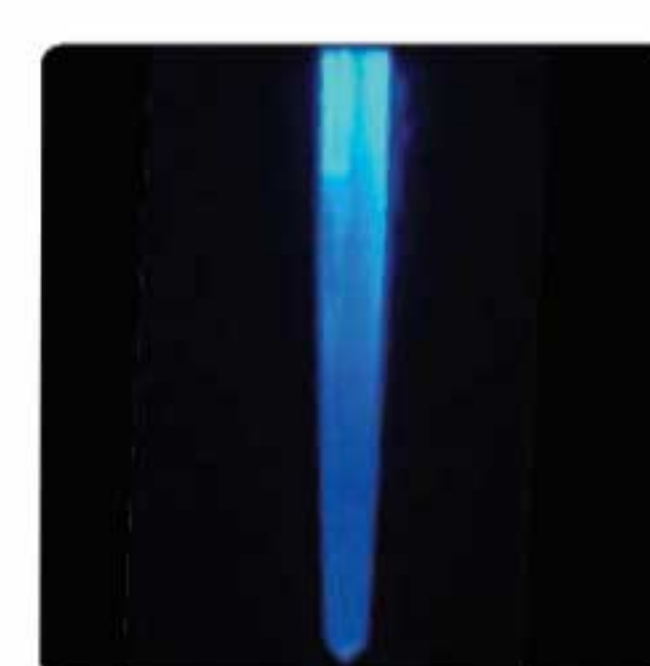
Light Transmission

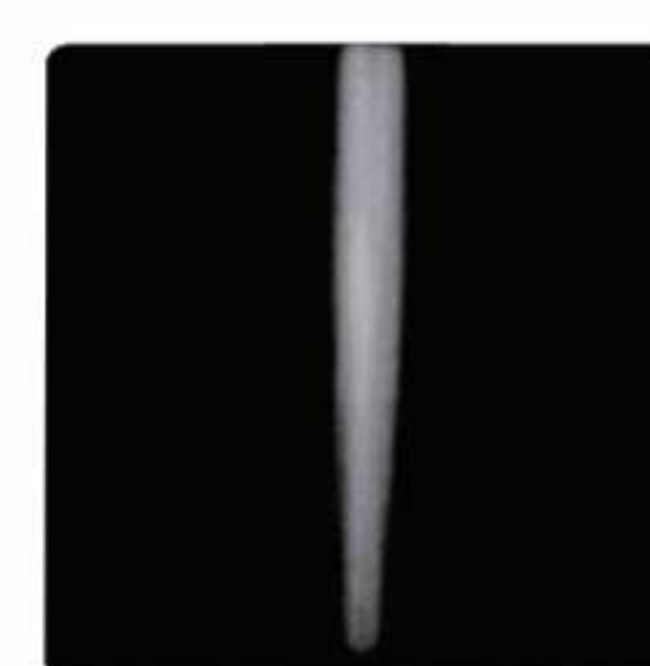Radiopacity

Scanning Electron Microscope

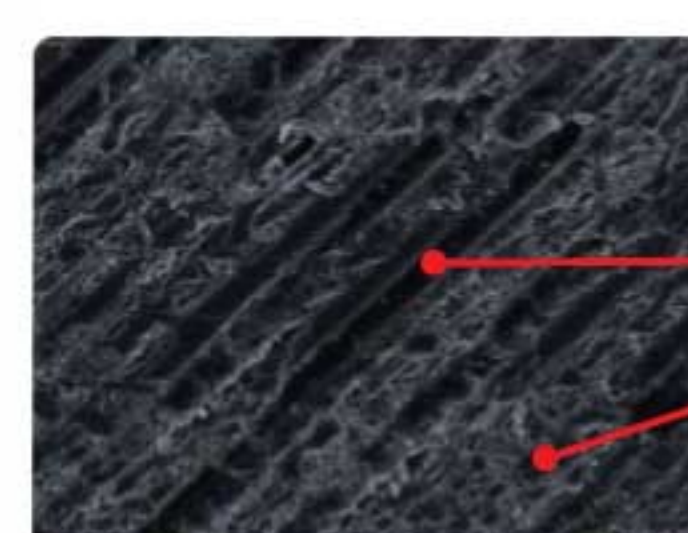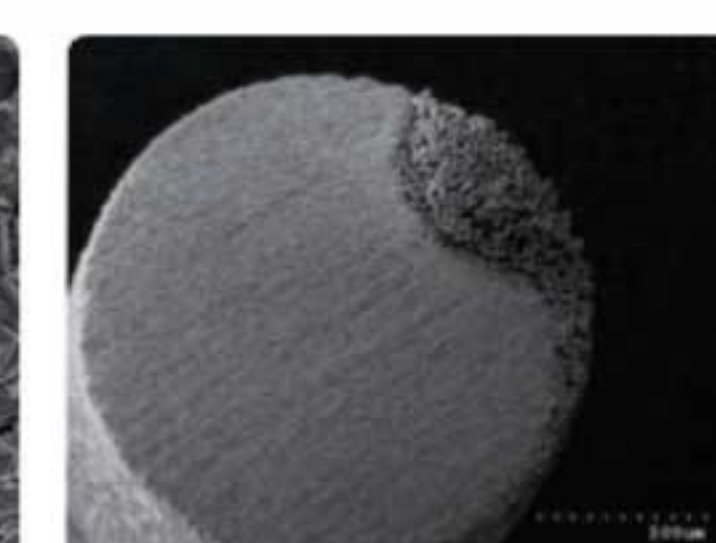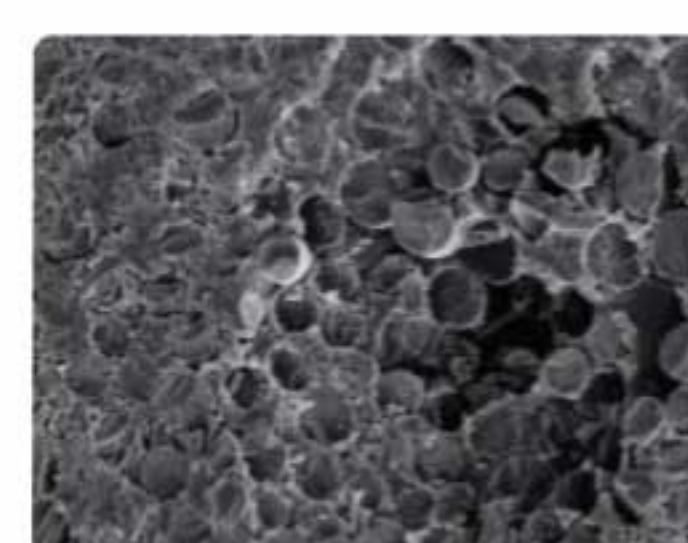Glass fiber
Composite resin matrix

Shelf life

- Three years from the production date.

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Applications

- Periodontal defects
- Intrabony defects
- Ridge augmentation
- Extraction sites(Implant preparation / placement)
- Sinus lift
- Cystic cavities

Advantages

- Excellent osteoconductivity
- Speed up the bone regeneration & heal activation
- Outstanding biocompatibility
(No biological rejection, No disease transmission)
- Accurate and problem-free implantation
- Easy handling
- Porosity: 55~60%

Implant Case


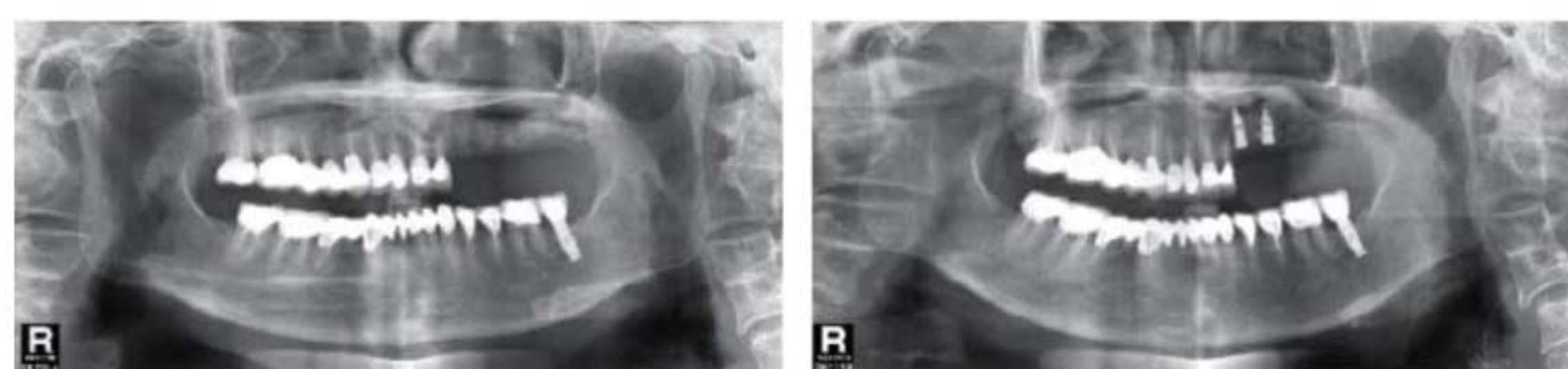
Fig. 1. Radiographic view of before implantation(Left), and after implantation(Right)

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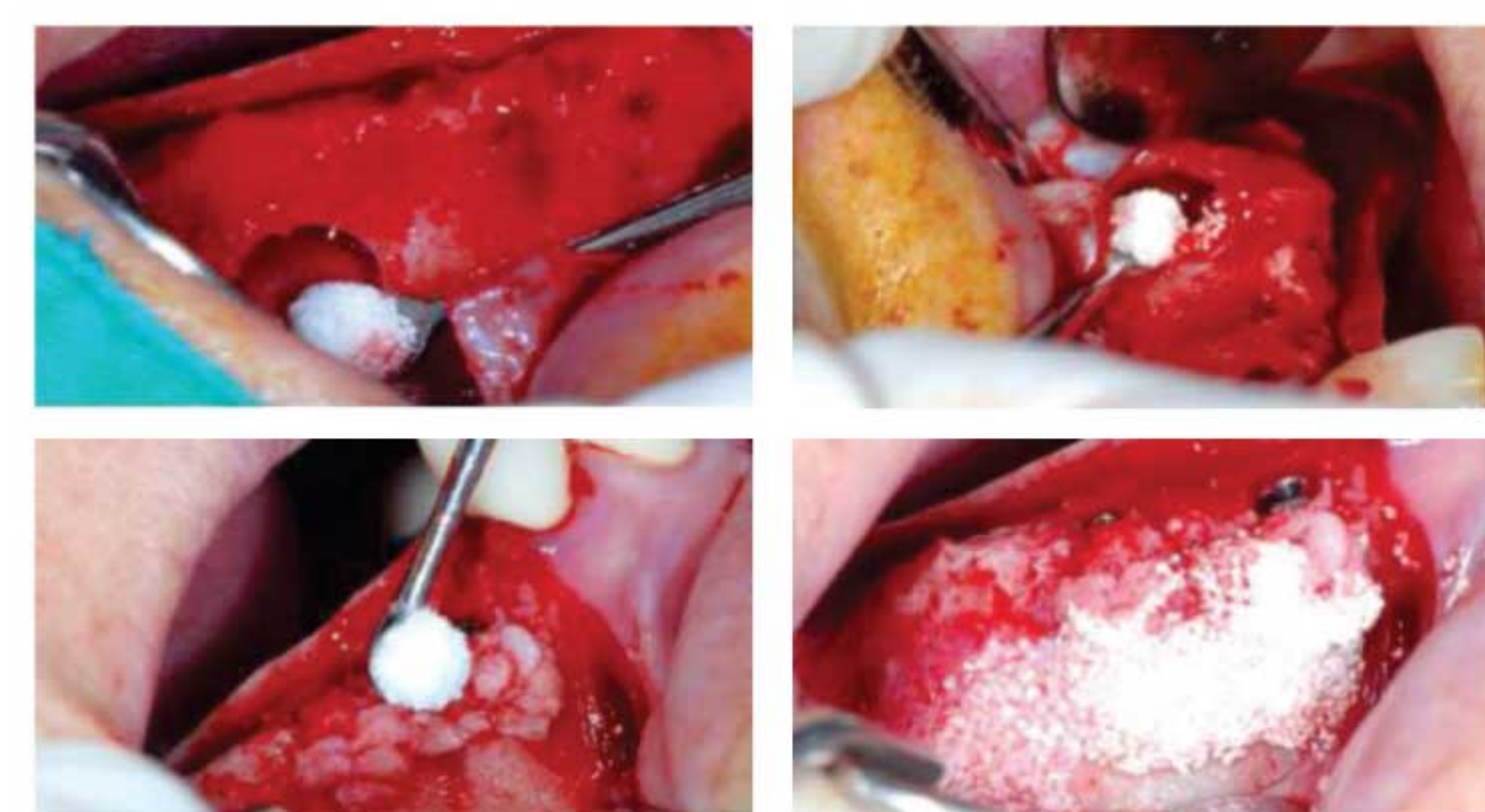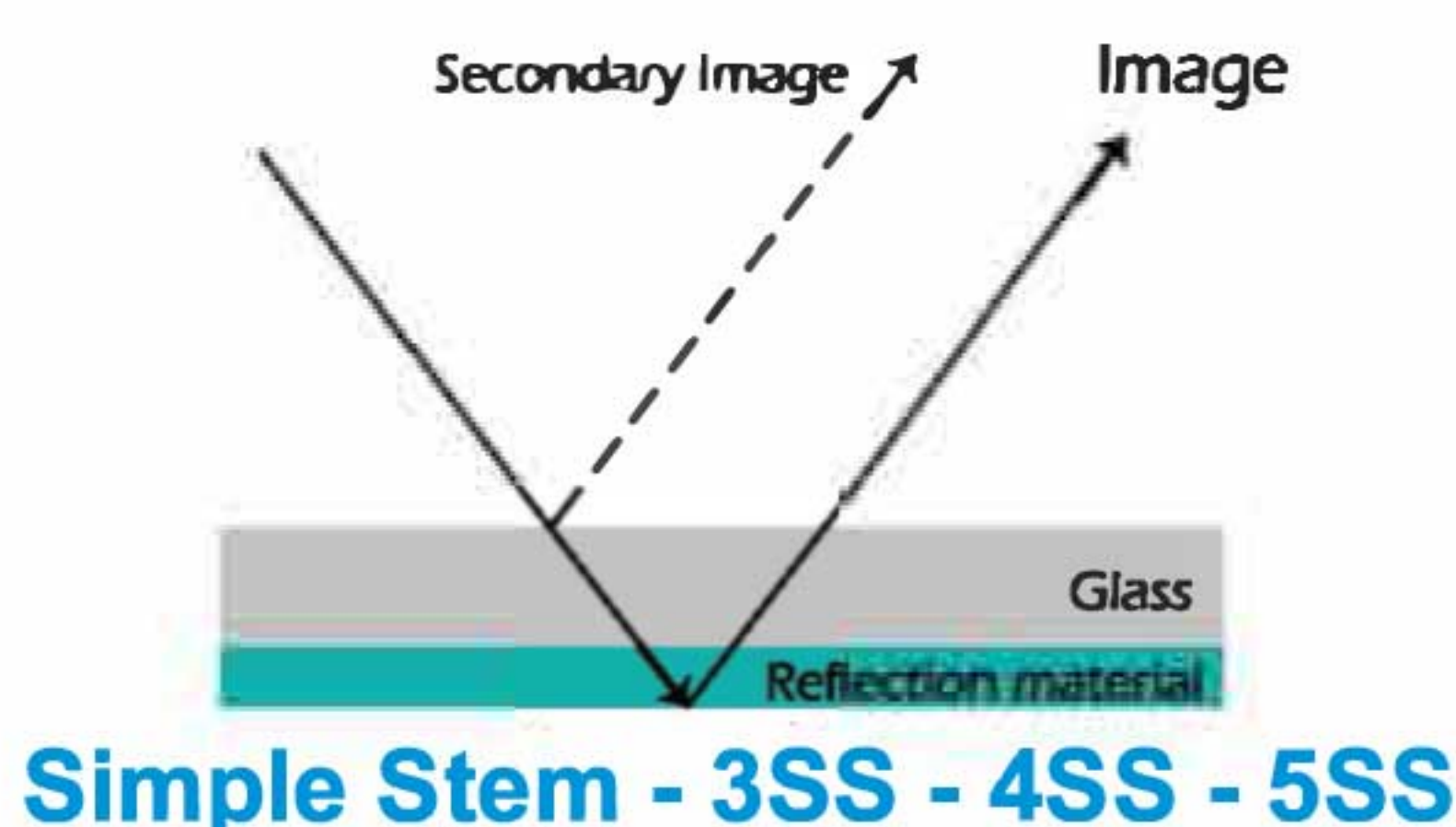
Advanced irrigation pump

The pump allows easy set-up of irrigation tubes and is extremely quiet during operation


Foot Control


- ① Coolant button(4 steps) ③ F/R change button
② Program change button ④ Foot pedal(speed control)

Optic

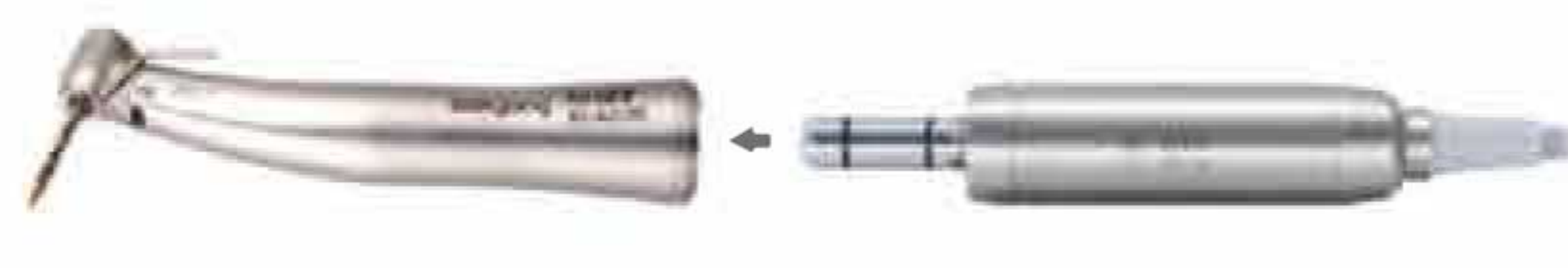Non-optic

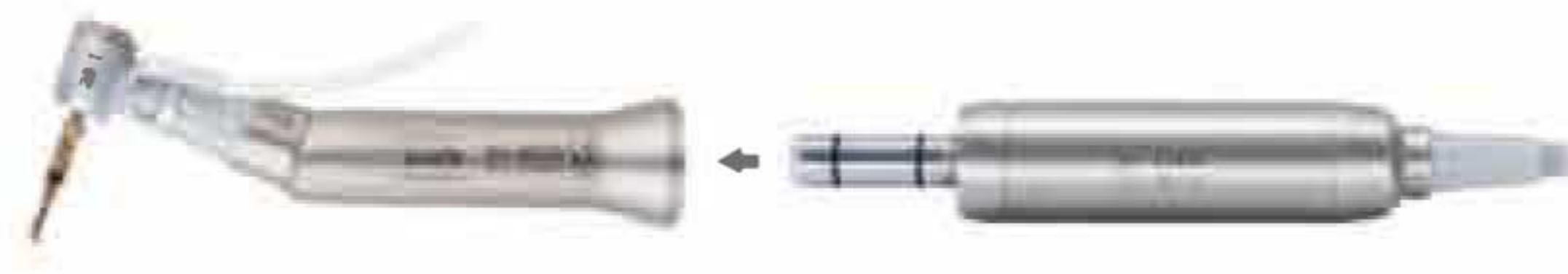MODEL Ki-20(Optic)

Specifications control box only

- Power supply : AC100-120 or AC220-240V 50/60 Hz
- Speed : 20-2,000 (20:1), 13-1,200 (32:1)
- Torque : 70Ncm (20:1), 90Ncm (32:1)
- Weight : 2.6 kg
- Dimensions : W240 x D225 x H108(mm)
- Contents : Control unit, Handpiece, Foot control, Handpiece stand

MODEL Ki-20(Non-optic)

Specifications control box only

- Power supply : AC100-120 or AC220-240V 50/60 Hz
- Speed : 20-2,000 (20:1), 13-1,200 (32:1)
- Torque : 70Ncm (20:1), 90Ncm (32:1)
- Weight : 2.6 kg
- Dimensions : W240 x D225 x H108(mm)
- Contents : Control unit, Handpiece, Foot control, Handpiece stand

KI-20 (Dual)

Implant Motor

**soeyang
krofit**
Made in Korea

Advanced irrigation pump

The pump allows easy set-up of irrigation tubes and is extremely quiet during operation.

Torque calibration

Calibrates the micromotor & the handpiece to the rotational resistance of each individual handpiece. During the operation, high torque accuracy offered.

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Data management

Treatment data can be easily accessed and downloaded using a USB stick. Files can be transferred and added to patient records.

*USB memory stick not included

MODEL New Ki-20 Advance (Dual)

Specifications control box only

- Power supply : AC100-120 or AC220-240V 50/60 Hz
- Speed : 10-2,000 (20:1), 7-1,200 (32:1)
- Torque : 70Ncm (20:1), 90Ncm (32:1)
- Weight : 3.893 kg
- Dimensions : W260 x D230 x H104(mm)
- Contents : Control unit, Handpiece, Foot control, Handpiece stand

i-ROOTS

Full Automatic Electronic Apex Locator

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i-ROOT is 5th generation Apex Locator by dual frequency type. i-Root can measure the apical constriction exactly without mismeasurement.

	Battery	Monier connection
i-Root S	Rechagable Battery	○
i-Root A	Rechagable Battery	×
I-Root 100	A2 Battery (3ea)	×

Advantages

- Best accuracy in any root canal condition (Dry, Bleeding, Wet, Saline, EDTA, NaOCl or Chlorhexidine etc...)
- Simultaneous visual and sound measurement
- Color display
- Stylish and ergonomic design
- Simple and easy to use
- Self test for the system checking
- Auto power off function
- Stable and durable in wide range of condition
- Computer monitoring & Patient Management

i-Root-S
(i-Root+PMS)

Specifications

Dimensions	W110 × D100 × H117mm
Screen	W79 × D70mm
Battery	Lithium-Ion, DC8.2V 2200mAh
Continuous use	60 hours
AC-Adapter	AC 100 ~ 240V, (50 ~ 60Hz)
Weight	390g

Components

- 1 i-ROOT main unit
- 1 Probe cord
- 1 File holder A
- 2 File holder B
- 5 Lip clips

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- 2.5D Patient Management Software, Fast, Reliable, Precise measurement results reduces need for X-Ray maximum portability for easy and quick working length gauging.
- Accuracy guaranteed by the latest multiple frequency technology meets hand-held size for the best efficiency of your endodontic treatment.

Advantages

- Accuracy in any root canal conditions (Dry, bleeding, wet, saline, EDTA, NaOCl, Chlorhexidine and etc)
- Color display
- Stylish and ergonomic design
- Auto power-off for saving battery life
- Stability and measurement reliability

Specifications

Rated voltage	DC 4.5 V (Alkaline battery size AAA DC1.5V x 3ea)
Input current	less than DC 60mA
Power consumption	less than 0.27 VA
Patient auxiliary current	less than AC 10uA
Oscillation frequency	500Hz±0.2Hz, 5kHz±0.002kHz
Protection type and level against electric shock	inner Power Type Equipment, Type B Equipment

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Endo E LED APEX

Endodontic Micromotor + Apex Locator

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Endo A LED APEX

Endodontic Micromotor + Apex locator

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Type per customer's request

Clinical Cases

Mandibular right 1st molar, mesial 4 orifices


Mandibular left 1st molar, 'S' curved distolingual canal

Mandibular left 2nd molar, C-shape canal

Downpack with Pack

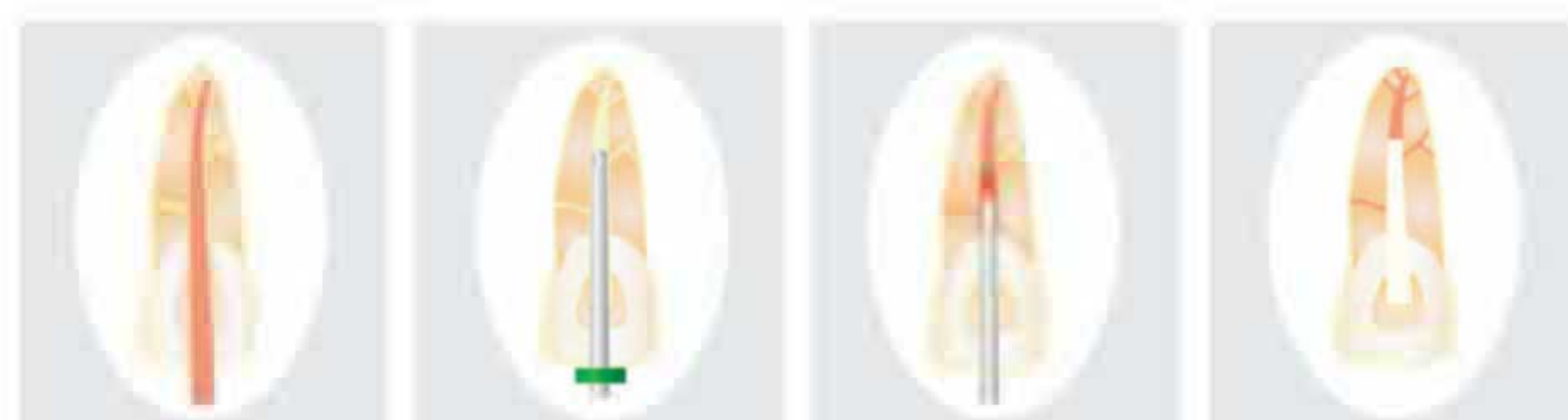1. Fit a master apical cone that matches the taper of the canal being obturated and check it makes a tug back.
2. Marked length select a Pack Tip and plugger 4-6mm shortly inserted than working length.
3. Cement master cone to 3mm, insert it to root canal and cut the excess over the orifice using the Pack Tip with gutta percha to 4-6mm short of the working length.
4. Harden the melted gutta percha with pre-selected plugger. Now apical area & accessory canals are perfectly obturated.

Backfill with Fill

5. After insert in the Fill Cartridge into the Root Canal, wait for 5 seconds after inserting Gun Needle to root canal.
6. Pull the trigger slowly and fill the gutta percha up to root canal. Fill Cartridge would be pushed by filled gutta percha naturally.
7. Harden the heated gutta percha with bigger plugger.
8. Now, the entire root canal is fully obturated.

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Powerful light intensity : 2,500mW/cm2

Effective cooling system with the anodized aluminum handpiece housing

Overheat protection sensor equipped

Mode selection : 1. Plasma mode(2,700mW/cm2) – 2 / 3 / 4 / Ortho sec curing

2. High Power mode(1,500mW/cm2) – 5 / 10 / 15 / 20 sec curing

Light intensity 8*8mm (2,500mW/Cm2)

Wavelength range: 400~490nm (peak : 405~465nm)

Dimensions Handpiece(155L*30W*23H)

Charger (φ120*85H)

Weight: Handpiece 120g (with battery w/o light guide)

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- Sterilize in autoclave

Repeat autoclaveable Guide-cap to prevent cross-infection and make curing efficiency safety.

- LCD Screen

LCD Screen indicates condition and state of device more easily to user.

- 330° head twist

330° head twist makes easy to cure anywhere inside of mouth.

- Interchangeable Battery

Easy-to Change Battery

Enough Battery capacity

Rechargeable Battery alone


- Beep Sound

Beep Sound make user's directive judgment to prevent misuse.

- Various curing mode

Choose light intensity on 4 kinds of modes at your convenience.

"selectable mode and time what you want."

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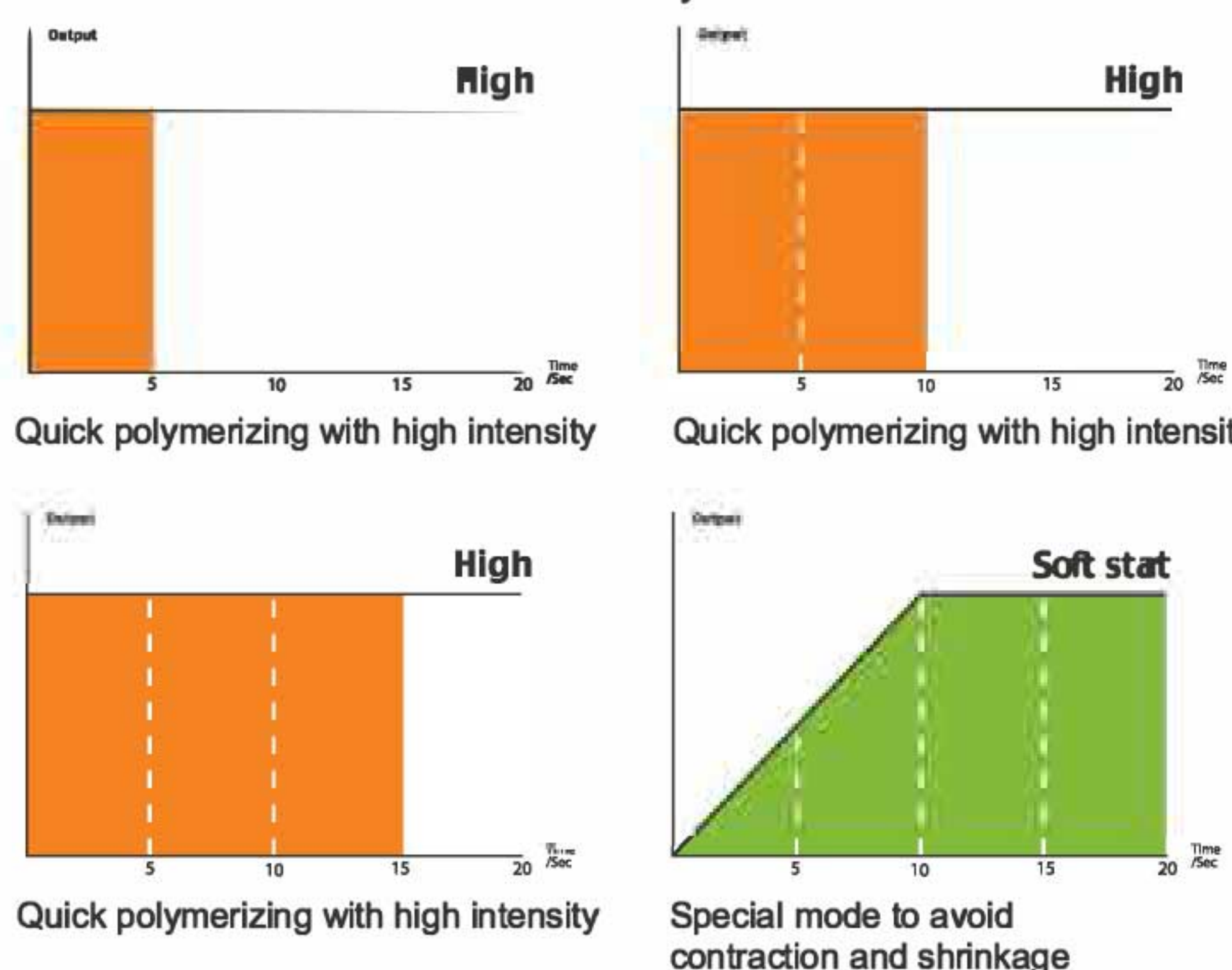E Driver

OSSTEM[®]
IMPLANT

Prevent orthodontic screw fractures, secure precise insertion path, minimize screw loosening and increase the convenience of abutment fastening

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All Prosthesis-related Screw Fastenings in One!

- Minimize fracture/loosening with precise fastening torque
- Easy access to posterior area by contra-angle structure
- Shorten chair time due to Rapid fastening

Torque Setting Guide

Prosthesis-related Screw	Torque (Ncm)	
	Mini	Reg
Cover Screw, Healing, Imp Coping	5	5
One piece ABT (Rigid, Solid ABT)	30	30
Two piece ABT (Transfer, Angled ABT)	20	30
Temporary Abutment	20	20
Cylinder Screw	20	20

From Orthodontic Screw Insertion to Removal!

- 15~35Nm torque allows insertion and removal, even at hard bone
- Torque setting function minimizes fracture/detachment of orthodontic screws

Orthodontic Screw	Torque (Ncm)	Speed (RPM)
Ø 1.4mm	15	
Ø 1.6mm	20	20 ~ 30
Ø 1.8mm	30	

K2 MOBILETM

Laser Engineering

SINCE 2002
HULASER
LASER SCIENCE FOR HUMANS

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
Application

SINCE 2002
HULASER
LASER SCIENCE FOR HUMANS

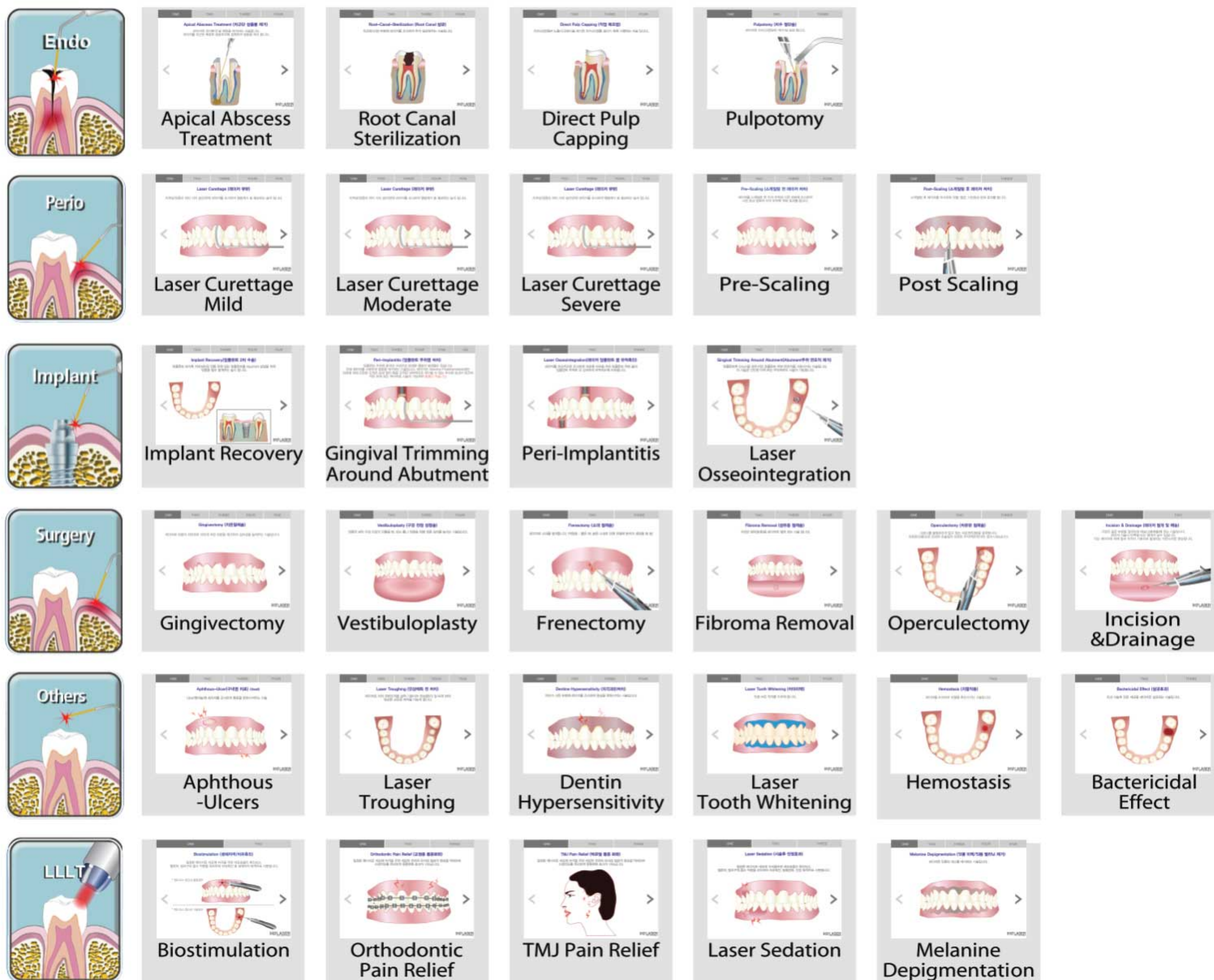Features

- Stable output power by HECO^{SA}™ Technology 3.5W Continuous wave 6W Pulse mode

- Intra-oral illumination

- No powercord and no footpedal
- Soft button technology
- Auto emission button


- Bright & clean OLED display
- Gravity sensor for both hands

- Ideal weight and size for handheld device -19cm, 108g

- Long lasting battery life -1700mAh
- Lithium-ion

- Various reusable tips (auto clvable) for endo, perio, surgery Tx, etc
- Least consumables
- Adjustable fiber length

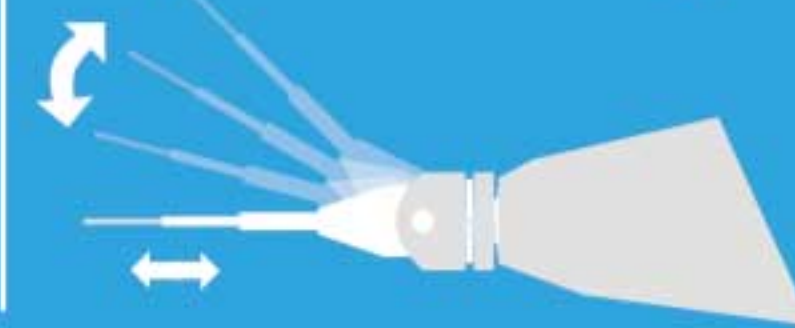- Head cap detachable and auto clvable

- World first movable fiber for least consumables and versatile usages

- Easy menu navigation

- Voice confirmation
- User interactive

- Fatigue-free and Ergonomic design (2014/2015 Good Design Award Winning)

- Fiber cleaver inside Cleaver blade to cut optic fiber for least consumables

Torque Handle

Machine Screw Drive Handle Attached to The Low-Speed Handpiece

OSSTEM[®]
IMPLANT

Tools

- Machine screw driver
- Corrective driver
- Bone profiler

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SlowJec


Minimize patient pain with automatic pressure adjustment function, excellent usability with one-touch operation

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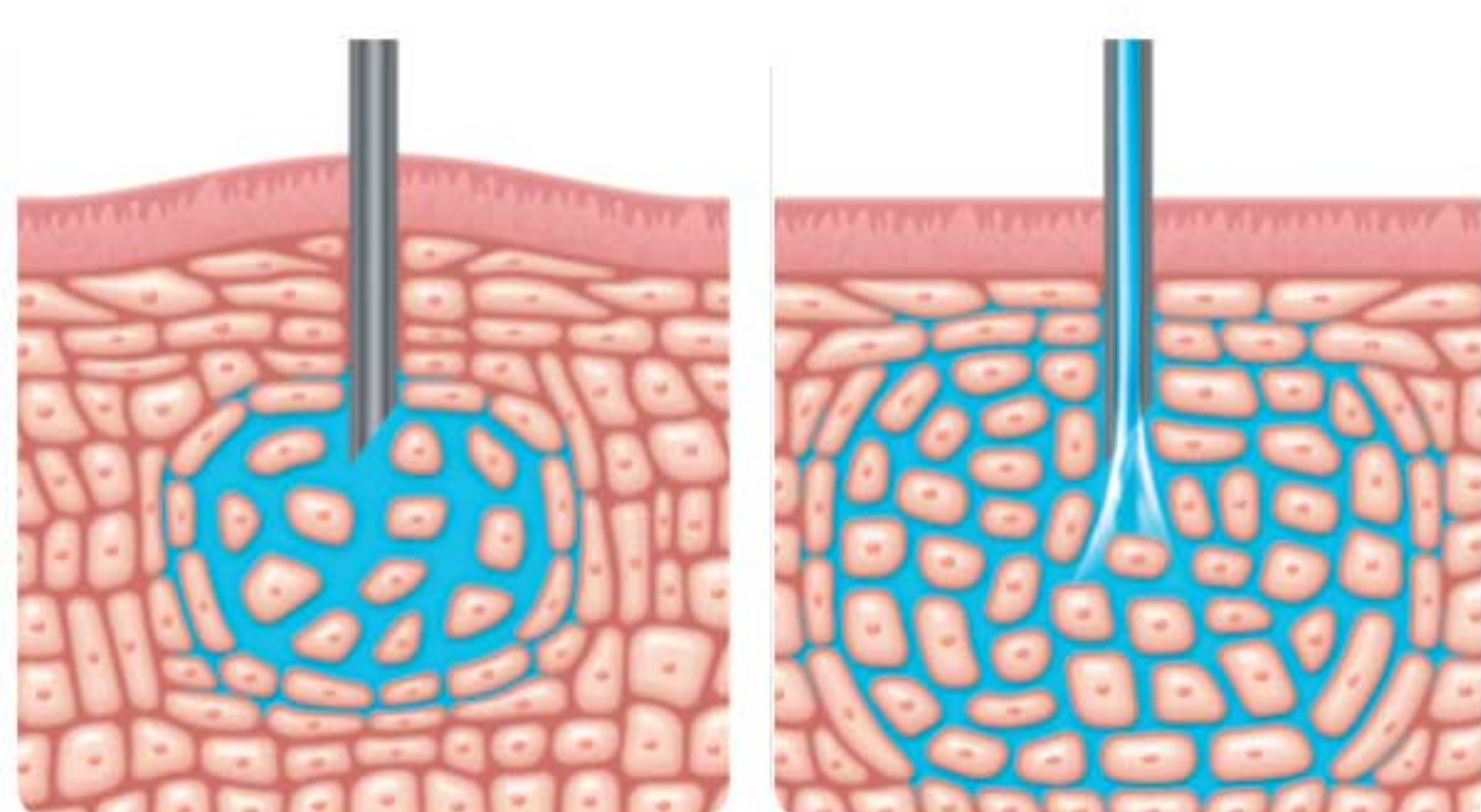Pause when over pressure is detected

Automatic injection when the pressure is decreased

Normal Syringe vs SlowJec

	Normal syringe	SlowJec
Injection time	With in 10seconds	1min 4seconds
Pressure	High	Low
Pain	Strong	Very weak

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Technical Specifications

Image sensor	CMOS (OV9712 or AR0330)	
Resolution	1.3M full HD (1280x720)	
Lens type	Liquid lens module	
Light source	White LED	
Focus range	5mm ~ 5cm	
Angle of view	70°	
Interface	Wireless	UWB
	Wired	USB 2.0
Operating voltage	Charger	DC 5V
	Battery (Wireless)	3.7V (li-ion)
	Receiver (Wired)	DC 5V

Dr.Cho's instrument Kit

OSSTEM[®]
IMPLANT

Optimal Implant Surgery Instrument Kit with Long-term Clinical Know-how

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1. Periosteal Elevator

- Comparing to Molt 9, with narrow width, fine and precise surgery possible
- The instrument lifting the mucosal periosteum after cutting the gingival tissue

2. Minesota Retractor

- Improving gripping area for long time use, which reduces hand fatigue (patent)
- The device that can secure the visibility by pulling the mouth, cheeks, etc.

3. Periosteal Elevator Selden

- The device that lifts and fixes the mucosal periosteum after cutting the gingival tissue

4. Needle Holder, Crile-Wood, TC

- Preventing the situation of being caught intermittently by bevel processing in the middle of the needle holder
- Tungsten carbide treatment in the Beak area enables retaining fine sutures, thereby increasing the convenience of the procedure.

5. Ochsenbein Chisel

- Useful when collecting a small amount of autogenous bone
- The blade is toward the front

6. Tissue Forcep, ADSON

- The instrument catching the gingival without hurting the soft tissue due to no protrusion on the inner surface.

7. Titanium Suction Tip

- High durability with titanium tip (inner diameter 3mm)

8. Dr.Cho's Instrument Pouch

- Using for the storage and sterilization of instruments
- Size 470X400mm

9. Extension Hose

- Hose for extension, connecting to unit chair suction
- Possible for autoclave sterilization (silicone material)
- Cleaning in transparent, Easy in dealing with

10. Extension Hose Adapter

- Adaptor to connect to unit chair suction

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

• Maxilla, Mandible single arch, Both arch, TMJ and etc.


FOV 8 X 8

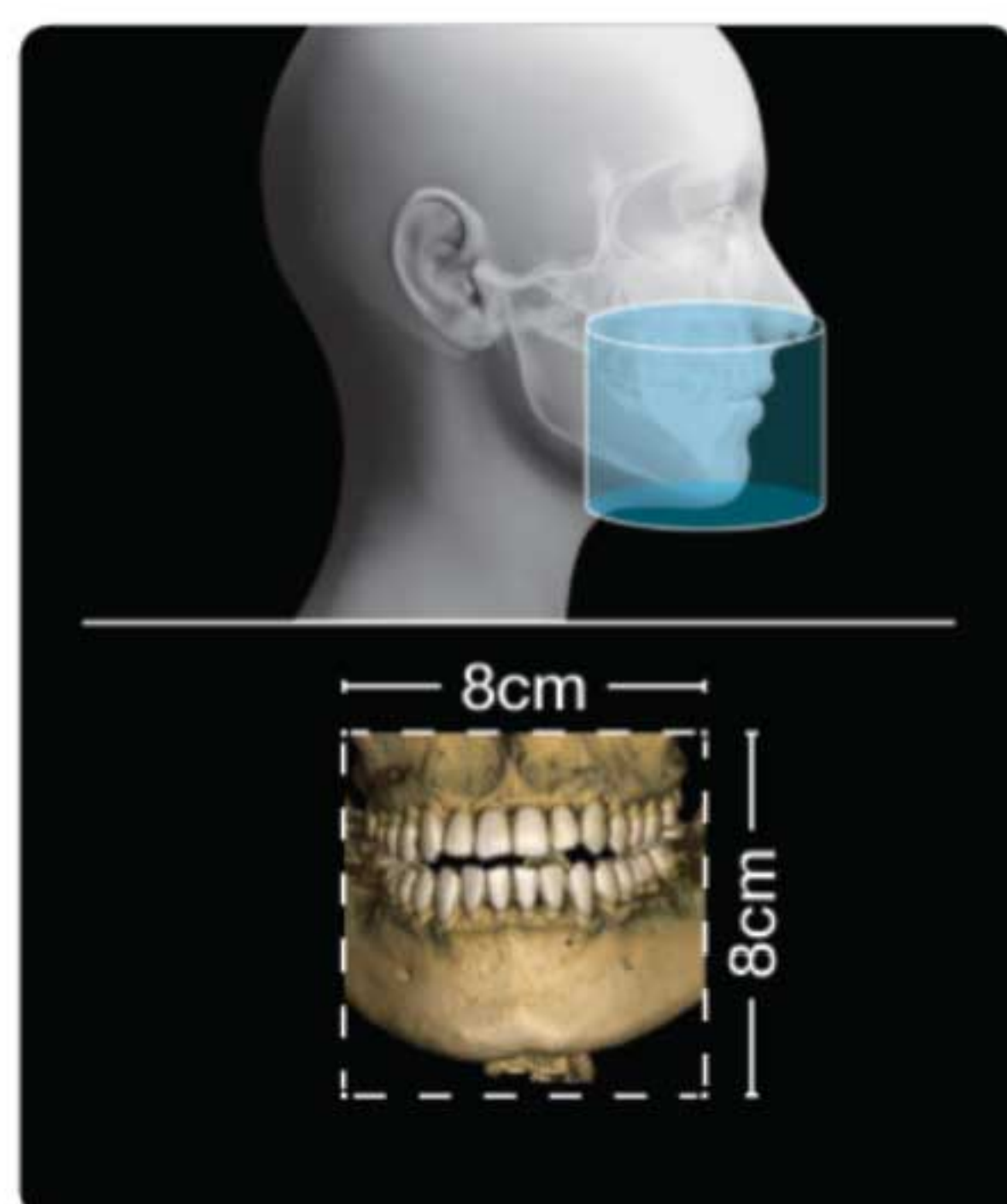
FOV 10 X 9

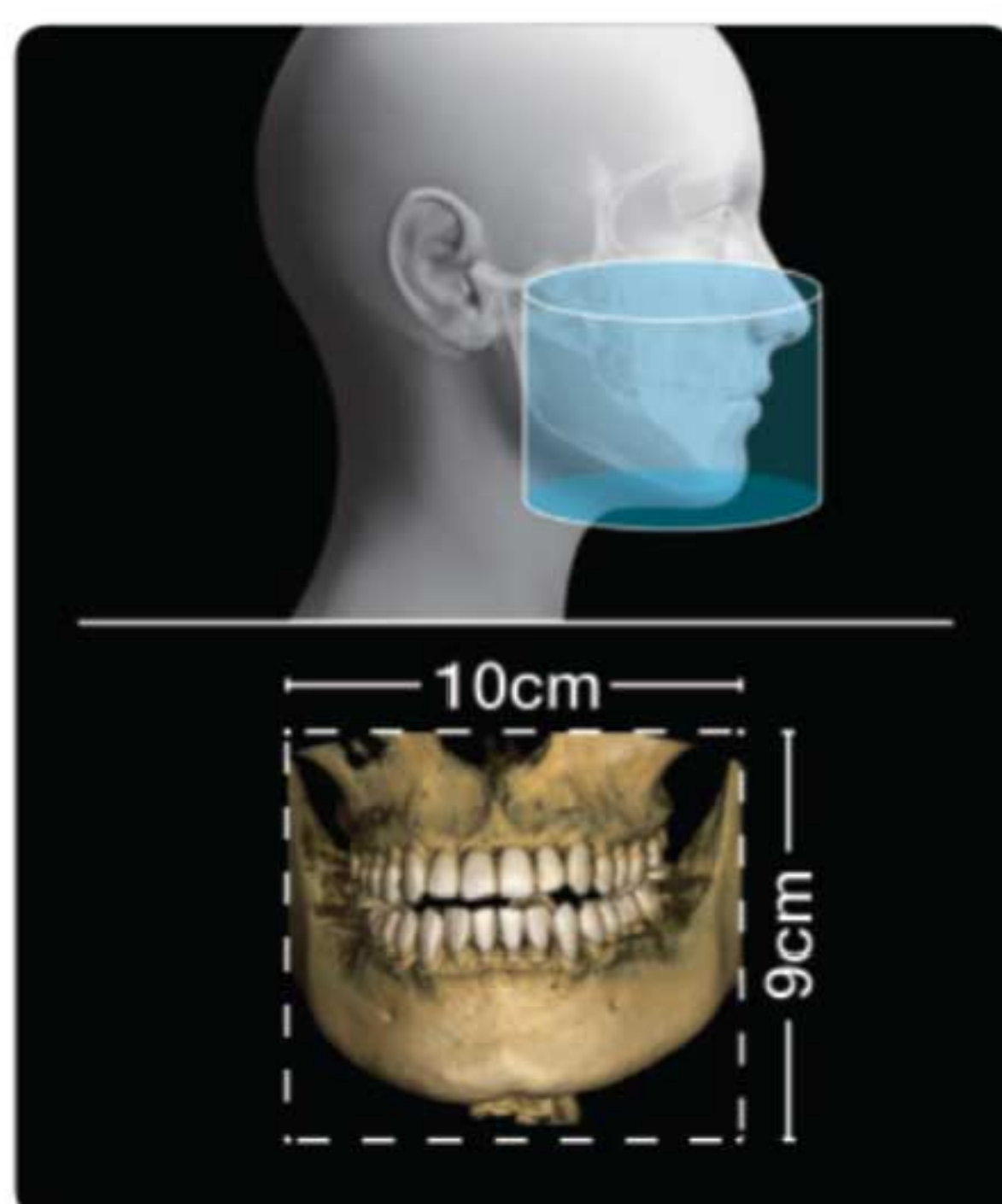FOV 12 X 9

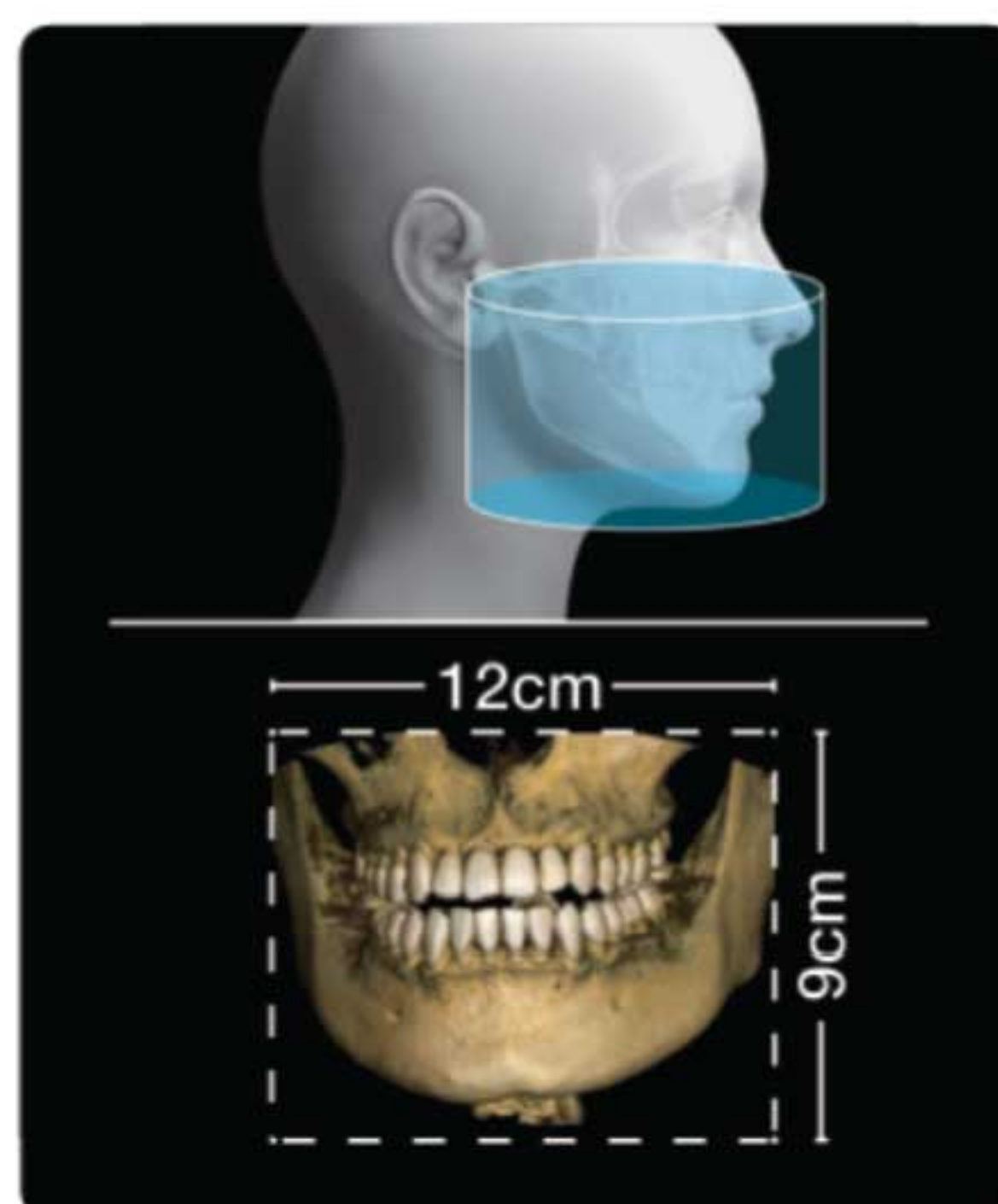FOV 15 X 9

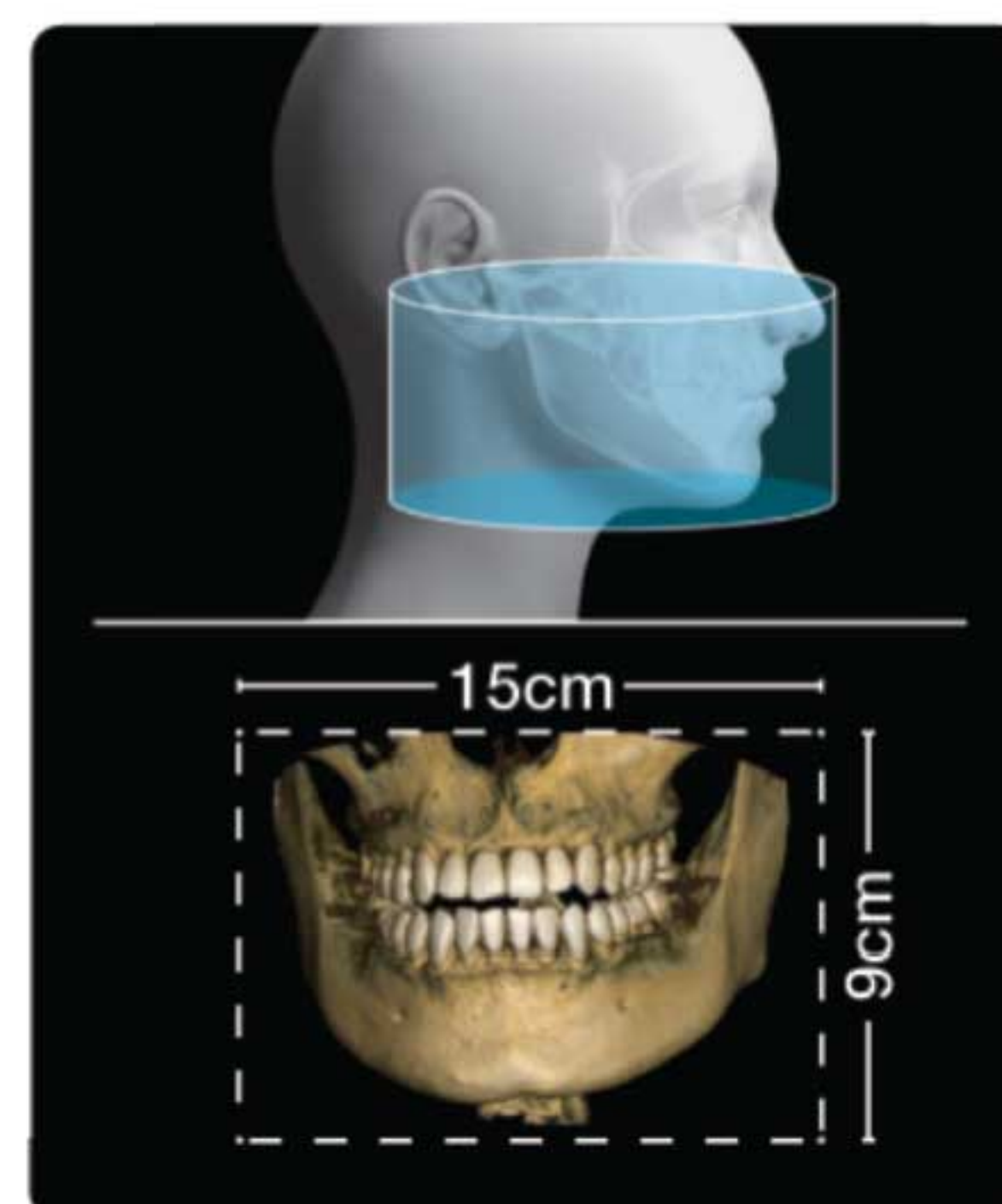Patient alignment laser

Mood light displaying status

Highly permeable support

Voice support

P2


Excellent Image Quality, Easy to Operate Panoramic X-ray

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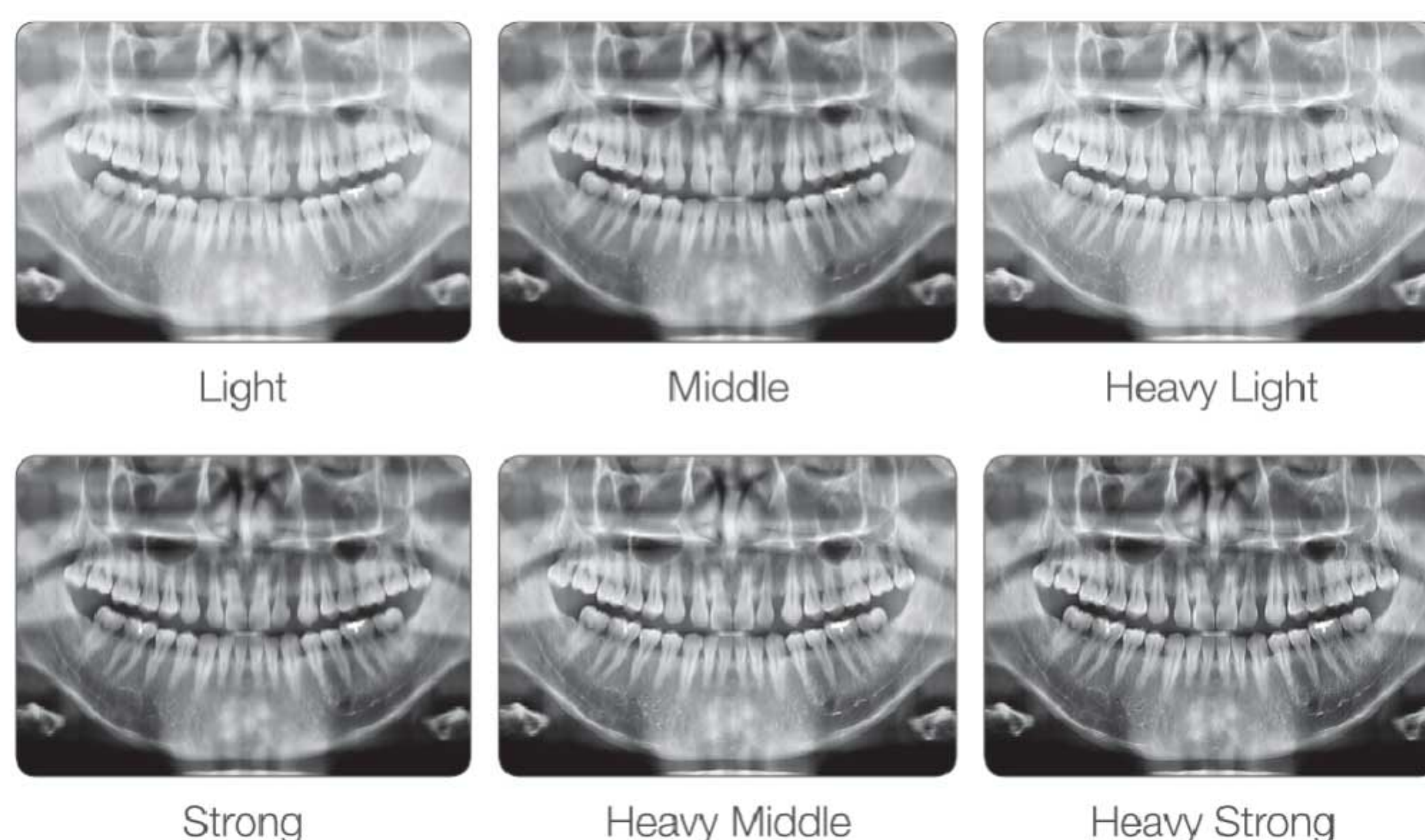R1

Sharp image, Easy-to-shoot standard X-ray

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

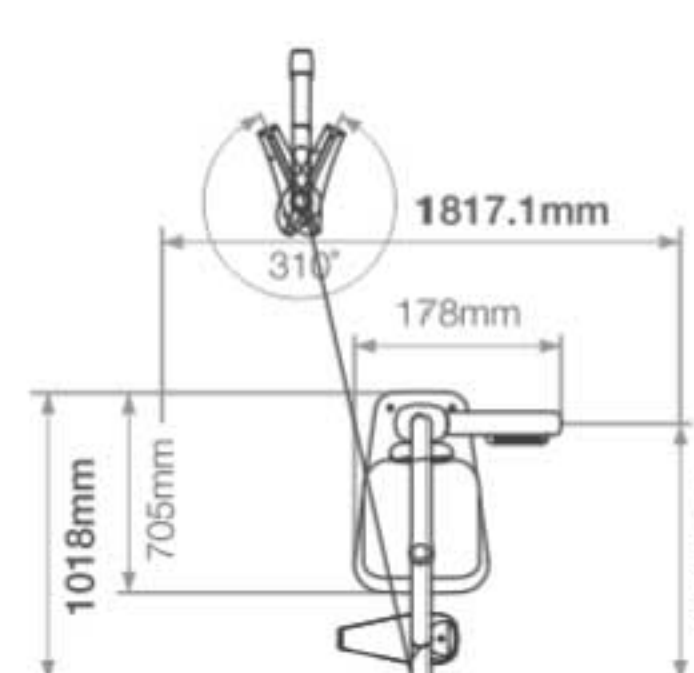
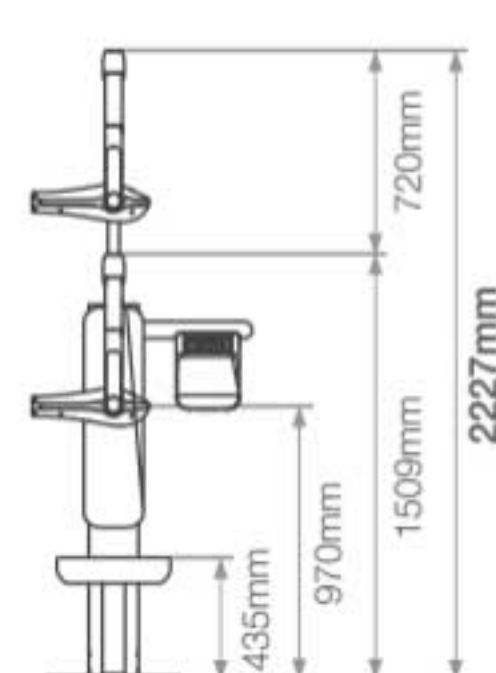
Wall Type

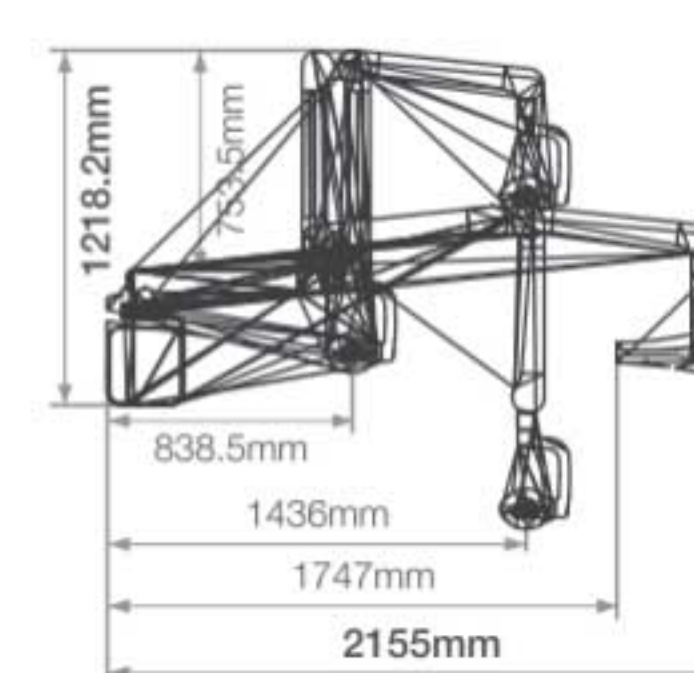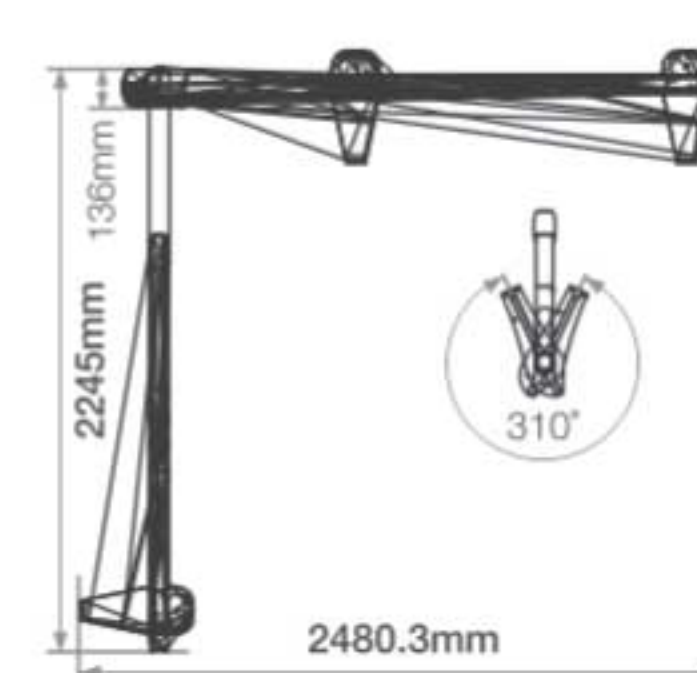BeauTis 15


Stronger whitening effect, Less toothache symptoms

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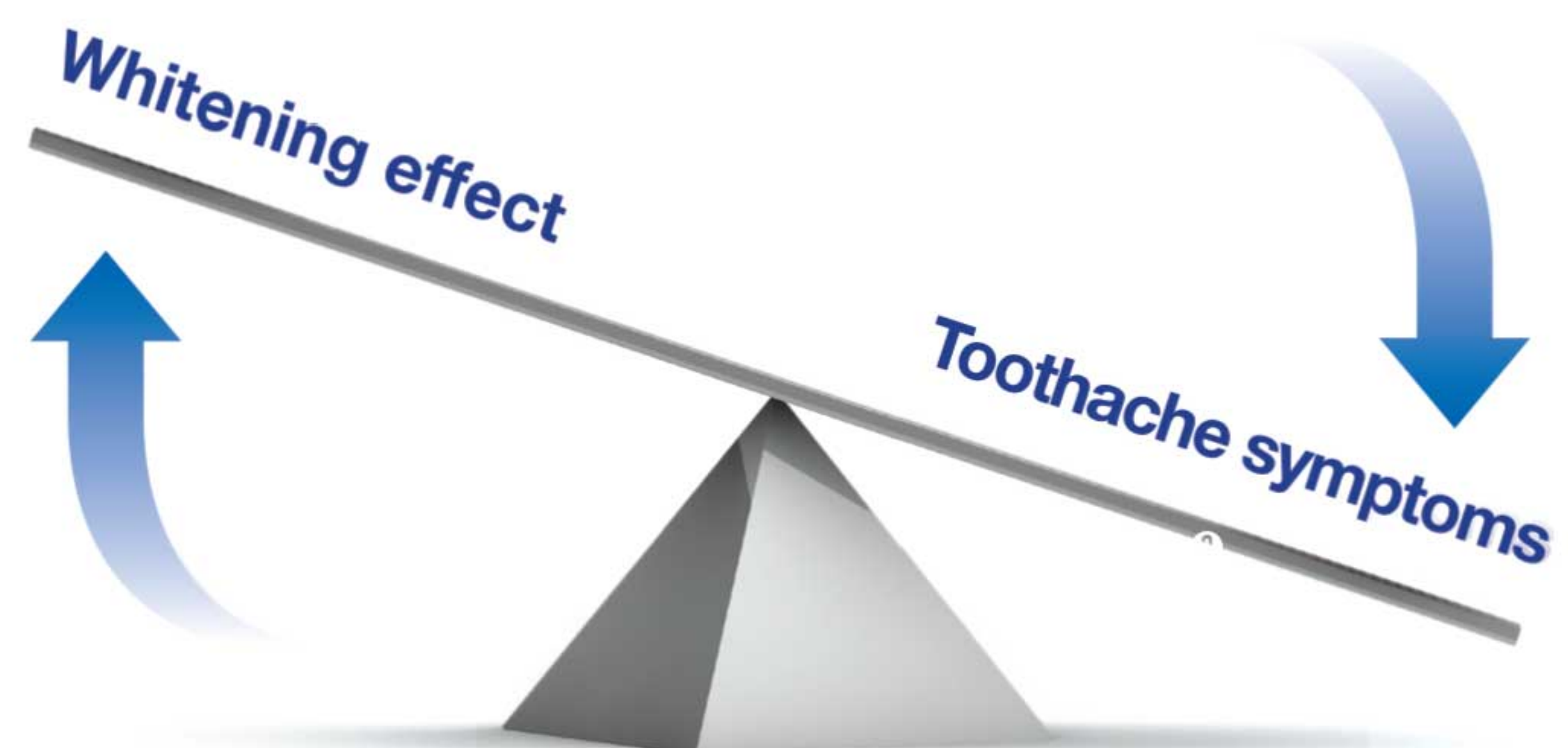Beautis White

Tooth whitening System that is Easy to Use in Daily Life

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Whitener Only (30minutes 2time a day, 2week use)

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Color Options

OSSTEM[®]
IMPLANT

K3 is available in a wide range of colorful options, so it blends harmoniously with any clinic interior.

1. Brown

2. Dark Blue

3. Red

4. Orange

5. Ivory

6. Light Green

7. Light Blue

8. Light Violet

The colors in the printed images may differ from those of the actual product

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Air Motor

OSSTEM[®]
IMPLANT

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Optic MODEL M500KL ORDER CODE P1005
Non-Optic MODEL M500K ORDER CODE P1007

Miniature Head : Power : 16 W, Speed : 380,000~450,000 min⁻¹, Head Size : ø10.6 x H 12.4 mm

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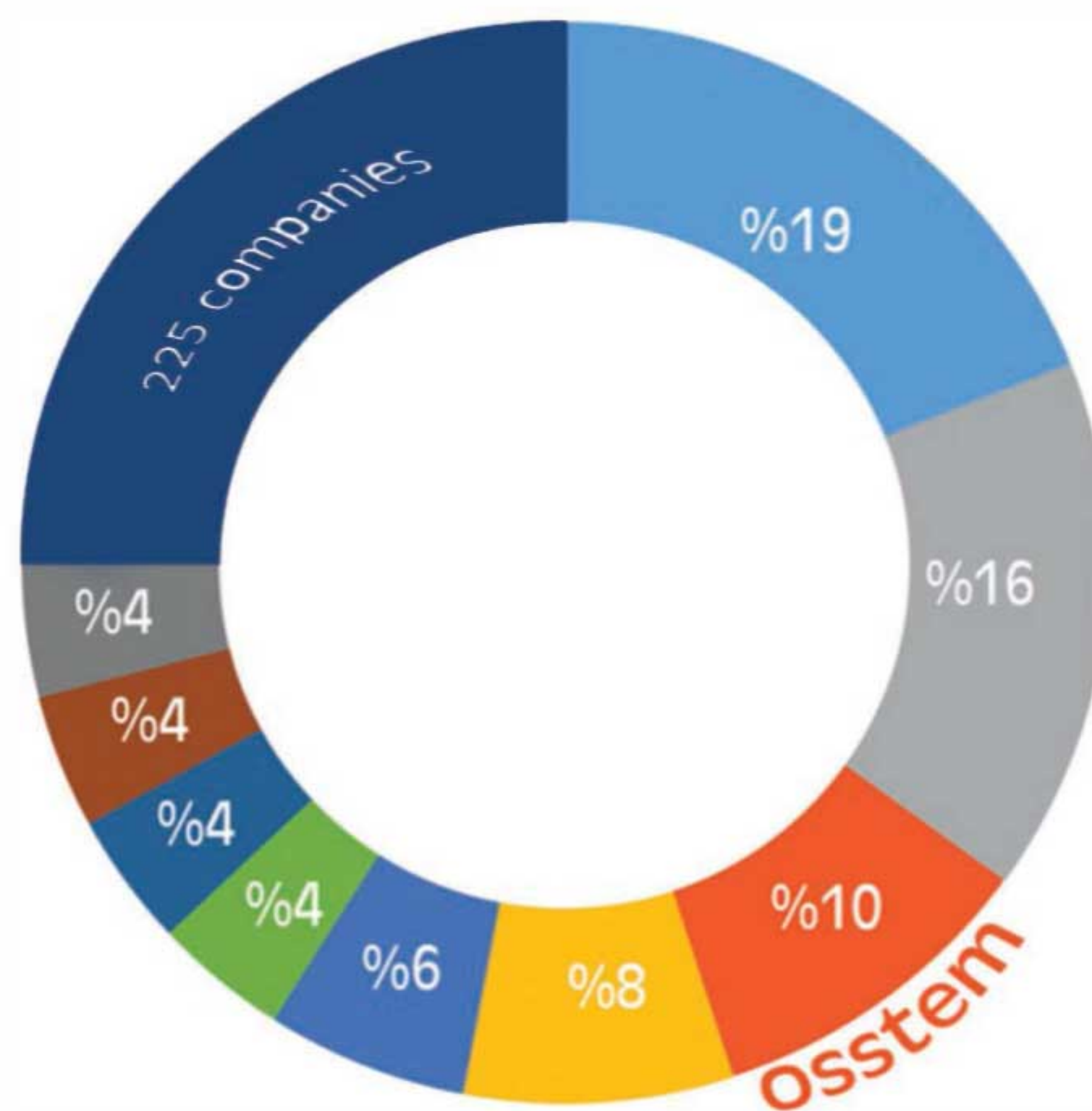
OSSTEM[®]
IMPLANT

Single Pitch Micro Thread

- 0.8pitch x 0.25depth x double lead
- Synchronized thread
- Optimized design for SA surface
- Reinforce fixture strength

Corkscrew Thread

- 0.8pitch x 0.5depth x double lead
- Powerful self threading
- Keep implant path
- Easy to change the path
- High initial stability
- Decrease sensitivity on drill size

Open Thread

- Prevents bone necrosis

Helix Cutting Edge

- Powerful self threading
- Easy path change

Apical design

- Self drilling ability
- Increased initial stability at fresh extraction socket


SA Surface

OSSTEM[®] SA Surface
IMPLANT

SA: Sand blasted with alumina Acid etched surface

- Optimal morphology: Combination of crater and micro-pit
- Optimal surface roughness : Ra 2.5 ~ 3.0
- Early cell response : 20% faster than RBM
- Early bone healing : 20% faster than EBM
- Early loading possible after 6 weeks of placement.

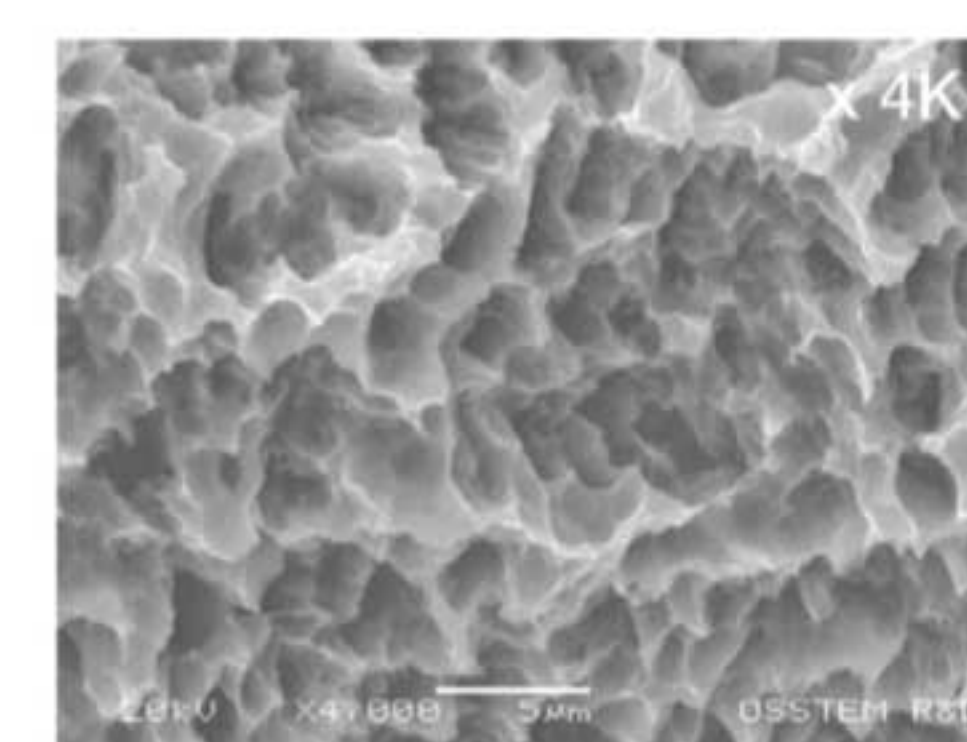SEM (Scanning Electron Microscope) image


CA Surface

OSSTEM[®] CA Surface
IMPLANT

Superhydrophilic SA surface encapsulated in calcium solution

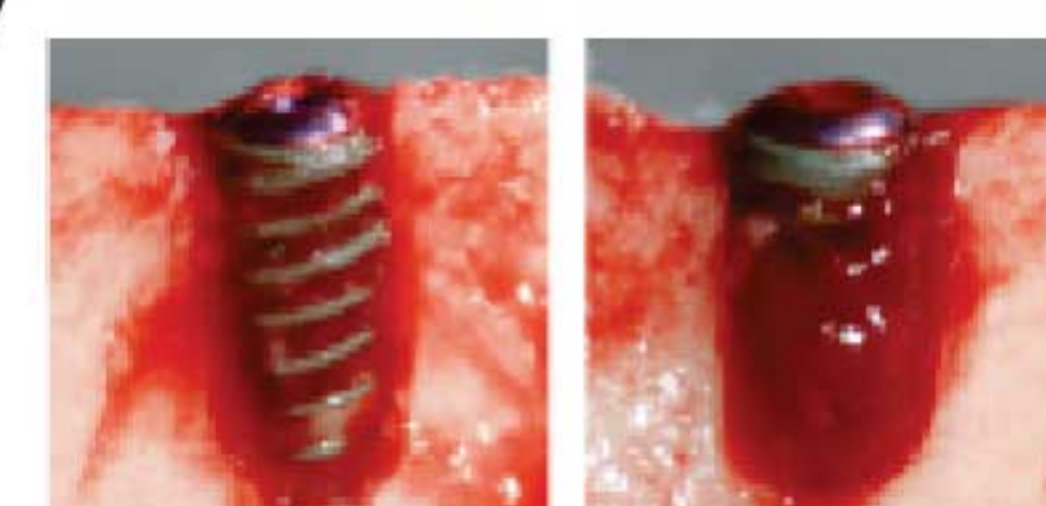- Maintains optimum surface identical to SA surface
- Surface activity maximized after encapsulated in calcium (CaC₂) Solution
- Increased ossification surface area through excellent blood wettability
- Improved bone reaction performance in the early osseointegration stage compared to SA surface

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