


Order Information

Description				
Kit	Small Kit	Post 15EA, Drill 3EA (Ø1.05mm, Ø1.25mm, Ø1.5mm)		
	Full Kit	Post 25EA, Drill 5EA (Ø1.05mm, Ø1.25mm, Ø1.5mm, Ø1.8mm, Ø2.2mm)		
Refill	Nexpost Refill	#0.1	Ø1.05mm (.035)	5EA
		#0.5	Ø1.25mm (.04)	5EA
		#1	Ø1.5mm (.06)	5EA
		#2	Ø1.8mm (.08)	5EA
		#3	Ø2.2mm (.1)	5EA
	Nexpost Drill Refill	#0.1	Ø1.05mm (.035)	1EA
		#0.5	Ø1.25mm (.04)	1EA
		#1	Ø1.5mm (.06)	1EA
		#2	Ø1.8mm (.08)	1EA
		#3	Ø2.2mm (.1)	1EA

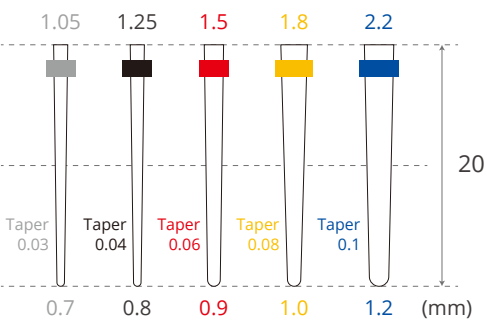

NexCore

Dual-curing core build-up material
Perfect Partner with Nexpost



Shade
Tooth, White, Blue

Components
Kit - 9g Syringe 2EA, Automix tip 20EA,
Intraoral tip 10EA, Endotip 10EA

Refill
- 9g Syringe 1EA

Meta Biomed Co., Ltd - Head Office
270 Osongsaengmyeong1-ro,
Osong-eup, Heungdeok-gu, Cheongju-si,
Chungbuk, Korea
Phone: +82 43 216 0433
info@meta-biomed.com
www.meta-biomed.com

Meta Biomed Inc.
3015 Advance Lane Colmar,
PA 18915,USA

Phone: +1 267 282 5893
america@metabiomed-inc.com
www.metabiomedamericas.com


Meta Biomed Europe GmbH
Wiesenstr. 35, 45473 Mülheim an
der Ruhr, Germany

Phone: +49 208 309 91910
europe@metabiomed-inc.com
www.meta-europe.com

DO-064-0

NexPost

Glass-fiber reinforced composite post

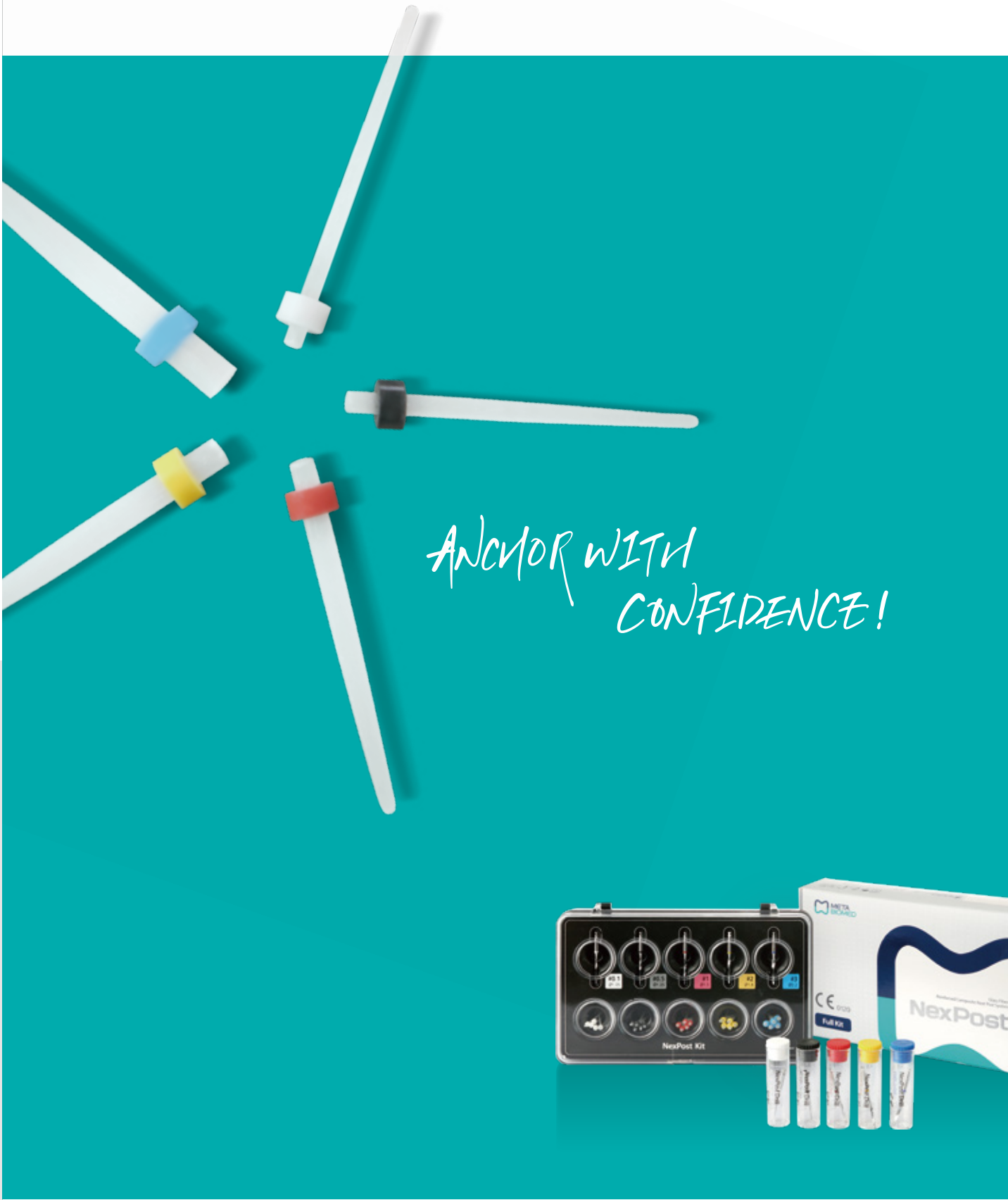

ANCHOR WITH
CONFIDENCE!



www.facebook.com/Metabiomed.Endodontics


Patented New-Technology

that can only be implemented by Meta Biomed

(Fiber post using Alkali-proof glass fiber and Making a method of it, No. 10-1313729)


High Strength

Nexpost is made of AR Glass fiber which has high strength, and it is cured in a vacuum state so that it can have higher flexural strength and bonding strength. This is our patented “New Technology” that can’t be imitated.

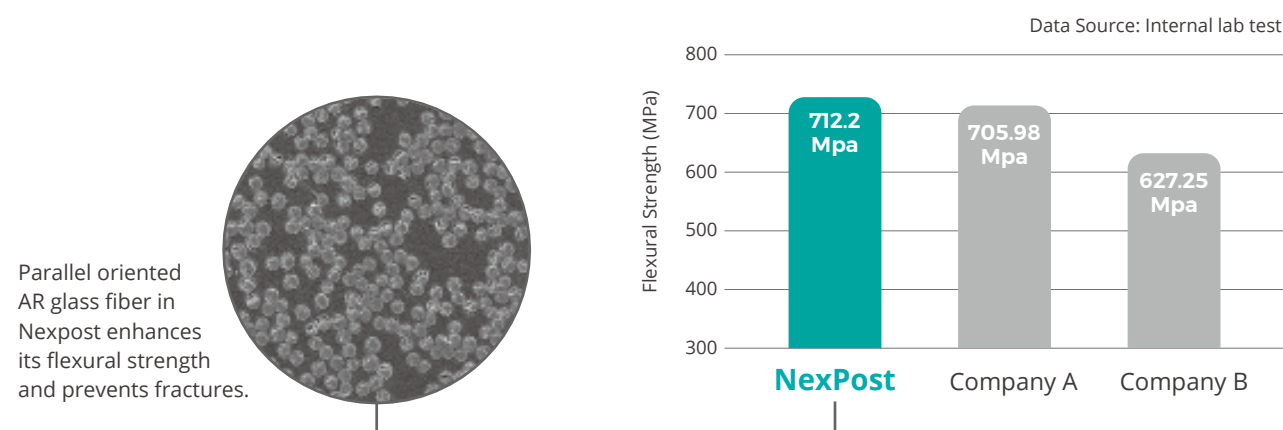

Similar Elasticity to Dentin


To prevent fracture of tooth and root, Nexpost has 21.7Gpa of flexural modulus (Dentin: 13.5-21.1Gpa) that can endure outside masticatory force (bite-force)

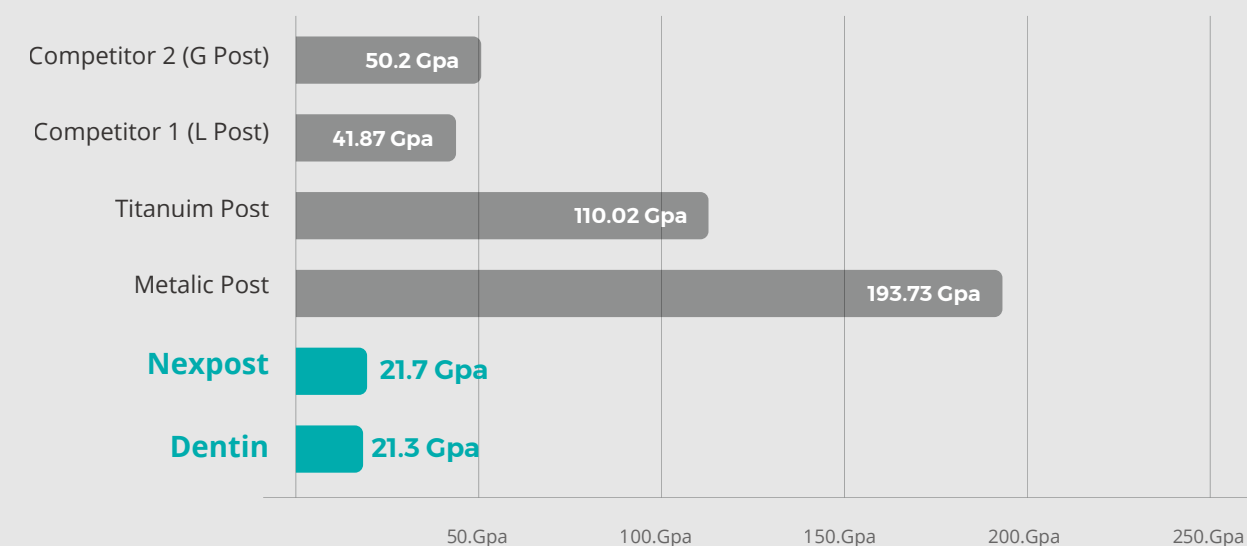

Data Source: Internal lap test & The flexural properties of endodontic post materials (Dominic A. Stewardson, Et al)


High Light-Transmittance

Thanks to our own patented “New Technology”, Nexpost has high light-transmittance allowing light-cure to be uniformly transmitted deep into the root canal.

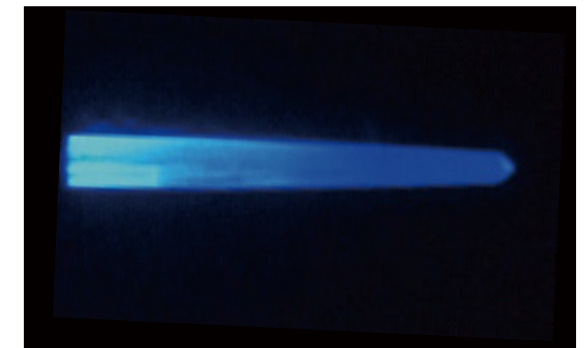


Excellent Radiopacity

Nexpost is made of Alkali-Resistant Glass fiber containing ZrO_2 (Zirconium Oxide) which has radiopaque properties, therefore Nexpost can have high radiopacity (390% Al)

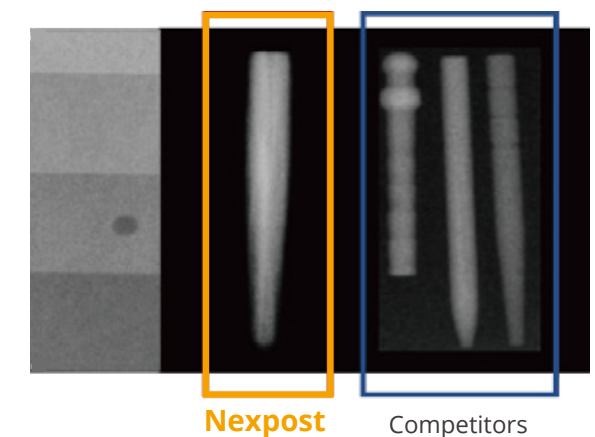

High Post-Core Bonding Strength

Especially with Nexcore

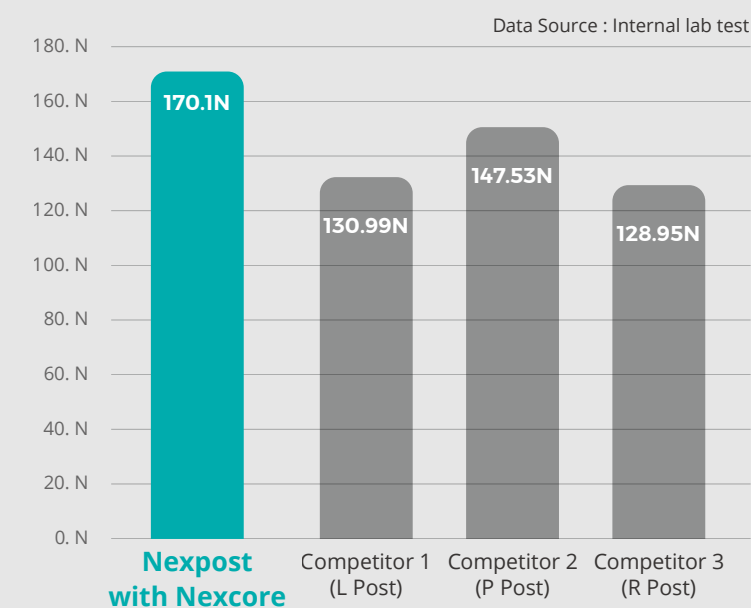

Surface of Nexpost is specially coated with Silane to enhance adhesion. Also Nexpost is made with Bis-GMA resin, which is also used in Nexcore (Competitors use epoxy resin) This process enhances adhesion to the core resin and causes exceptional compatibility to Nexcore.

